
Health o meter®

BODY FAT & HYDRATION MONITORING SCALE

User's Manual
Model BFM582

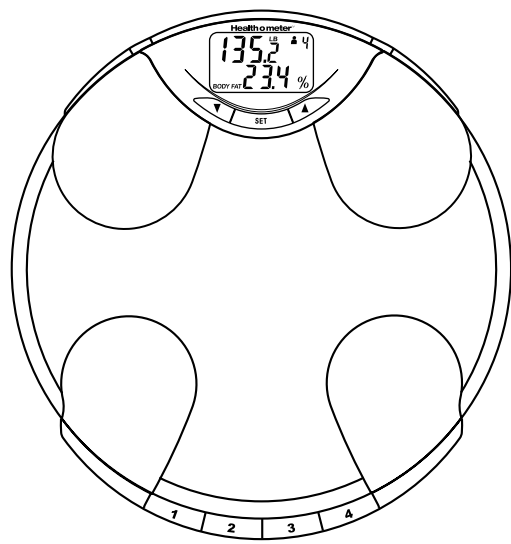PLEASE READ CAREFULLY AND SAVE

For information on service or warranty, contact Sunbeam Products, Inc. at 1-800-672-5625 or go to www.healthometer.com.

To assist us in serving you, please provide model number and date of purchase.

Monitoring your Weight, Body Fat Composition and Hydration Levels

- * Make sure you have programmed your personal data.
- * Place the scale on a flat, level surface. A hard surface is best to achieve maximum accuracy. Placing the scale on a rug or carpet may affect the scale's weighing accuracy.

1. Press your personal user number with your toe. Your personal profile will appear on the screen, along with last weight, body fat and hydration measurements.
2. Wait until "0.0" and "0.0%" appears on the display.
3. Stand on the scale (bare feet only!) so the balls of your feet are on the upper silver pads, heels are on the lower silver pads.

NOTE: For body fat and hydration estimates, always measure with bare feet.

4. During the weight measurement, the display will flash "-.-" and your weight will appear. Next, the display will flash "-.-", your body fat percentage and your hydration level will be shown.
5. The scale will turn off automatically after a few seconds.

NOTE: It is normal for your weight to vary during a day and from one day to another, slight changes in weight are generally normal.

Weight-only Usage

It is not necessary to have the user's personal data programmed into the scale to obtain an accurate weight-only reading.

- * Place the scale on a flat, level surface. A hard surface is best for maximum accuracy. Placing the scale on a rug or carpet may affect the scale's weighing accuracy.

1. Press on the scale platform with the ball of your foot to activate the scale.
2. Wait until "0.0" is displayed.

WELCOME

Congratulations on your purchase of this Health o meter® wellness monitoring product. It's so much more than just a scale. It helps monitor the weight, body fat composition and hydration level of the body for up to four different users. It's easy to program and use. Once your data is stored, your scale is ready to use. It's that simple!

ABOUT BODY FAT

The Health o meter® Body Fat Monitoring Scale can help you in achieving your physical best by providing three critical measurement estimates: BODY FAT COMPOSITION, BODY WEIGHT and BODY HYDRATION LEVEL.

Traditional bathroom scales, indicate your overall weight—both lean mass and body fat together. The Health o meter® Body Fat Monitoring Scale lets you track your weight, body fat and hydration separately.

IMPORTANT!

Please read all instructions to familiarize yourself with features and operations before programming the scale.

CAUTION!

- The use of this scale is not recommended for persons with an electric implant (heart pacemaker).
- The scale will not calculate body fat unless you have bare feet.
- For an accurate body fat estimate, do not bend your knees, and keep your legs/thighs apart and straight when standing on the scale.
- The body fat function is not recommended for use by pregnant women, persons with fever, swollen legs or other edemas, as well as a person who is over hydrated or dehydrated.
- Body fat percentages may be slightly higher or lower for children under 16, persons with diabetes and other health conditions.
- Women naturally store approximately 5% more body fat than men.

3. Step on the scale and stand still. The display will flash "-.-" and your weight will be shown.
4. The scale will turn off automatically after a few seconds.

NOTE: It is recommended to weigh yourself at the same time of the day. It is normal for your weight to vary during a day and from one day to another, slight changes in weight are generally normal.

Body Hydration

Body hydration percentage indicates the amount of fluids you have in your body. Water in your system is important because it regulates body temperature, converts food into energy and protects vital organs. Individuals with more body fat have proportionately less total body water and are more susceptible to fluid imbalances that cause dehydration.

When tracking your body hydration level you should consider the following:

1. Take your hydration measurements the same time of day.
2. Monitor your readings to establish a personal body hydration range.
3. Do not drink any caffeinated or alcoholic drinks before taking a measurement.
4. Drink water if dehydration is detected.

Extracted from Dialysis & Transplantation, © 2001; Sports Nutrition For Young Adults: Hydration; Crystal Springs: Hydration Information; The National Women's Health Information Center: 4woman.org

Please consult with your physician to determine the appropriate hydration level for your body.

CAUTION!

If "E" appears on the display, an error has occurred. Step off the scale and carefully repeat the steps again. The "E" will appear if:


1. You step on the scale before "0.0" displays.
2. The weight is greater than the maximum capacity (330Lb / 150kg). NEVER apply a load beyond the maximum capacity.
3. The scale cannot measure body fat.

OPERATING INSTRUCTIONS


Installing the Batteries

1. Locate the battery cover on the back of the scale.
2. Remove the battery cover and insert 4 AA batteries.
3. If "Lo" appears in the display, replace the batteries.

Note: Do not use rechargeable batteries.

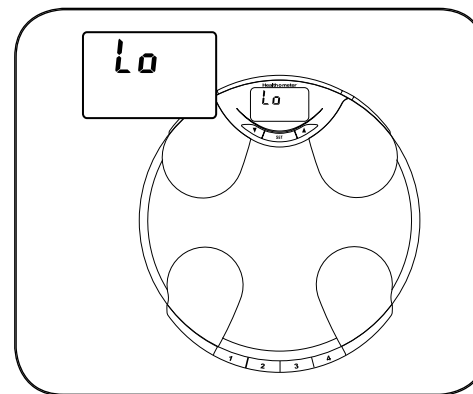
Keys and Icons


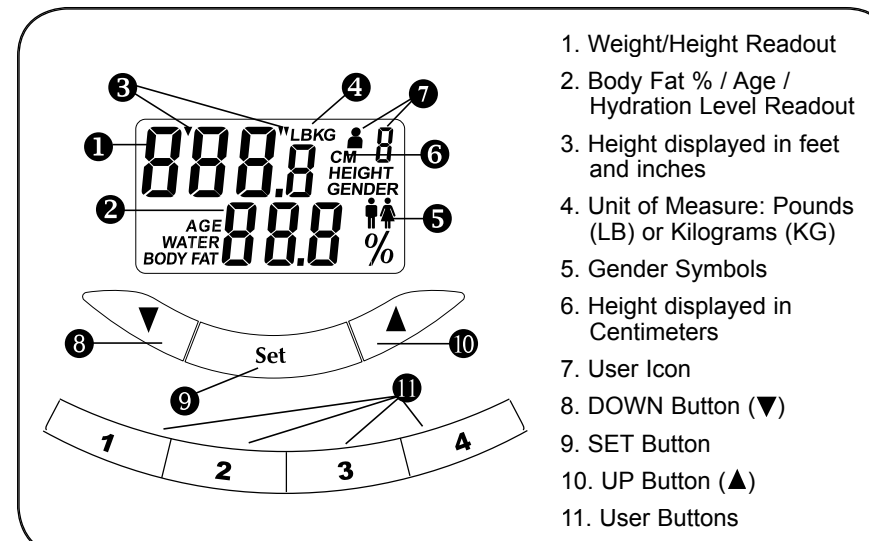
1. Weight/Height Readout
2. Body Fat % / Age / Hydration Level Readout
3. Height displayed in feet and inches
4. Unit of Measure: Pounds (LB) or Kilograms (KG)
5. Gender Symbols
6. Height displayed in Centimeters
7. User Icon
8. DOWN Button (▼)
9. SET Button
10. UP Button (▲)
11. User Buttons

Body Fat Composition Ranges

Men's Body Fat Composition Ranges

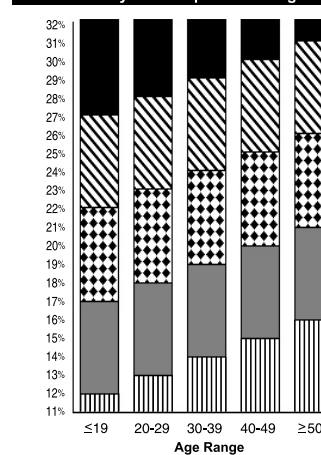Women's Body Fat Composition Ranges

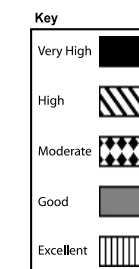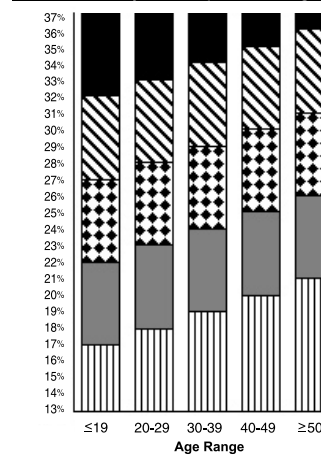Extracted from Principals & Labs for Fitness and Wellness, 1st Edition ©1999. To determine the percentage of body fat that is appropriate for your body, consult your physician.

Caring for your Scale

Periodically wipe the scale surface down with a DRY cloth to remove dust.


CAUTION! Weight exceeding the capacity (330 lb/150 kg) may damage your scale. This product is designed for personal use only. Not for commercial or industrial use. This unit is not waterproof; avoid contact with excessive moisture.

First Time Programming

On the back of the scale platform, use the selector switch to set the unit of measurement to "LB" (pounds) or "KG" (kilograms).

Turn the unit over and follow these steps to program your Body Fat Monitoring Scale:

1. Press and hold the desired User Button until the HEIGHT icon starts flashing.
2. Press the ▲ or the ▼ button to select the desired FEET (or CENTIMETERS in KG mode); then press the SET button.
3. Press the ▲ or the ▼ button to select the desired INCHES; then press the SET button.
4. Press the ▲ or the ▼ button until your correct AGE is shown on the display; then press the SET button.
5. Press the ▲ or the ▼ button to select your gender (MALE or FEMALE); then press the SET button.
6. Once you have entered your information, wait until the display shows "0.0". Your profile is now programmed into the selected user button and the scale is ready for use.

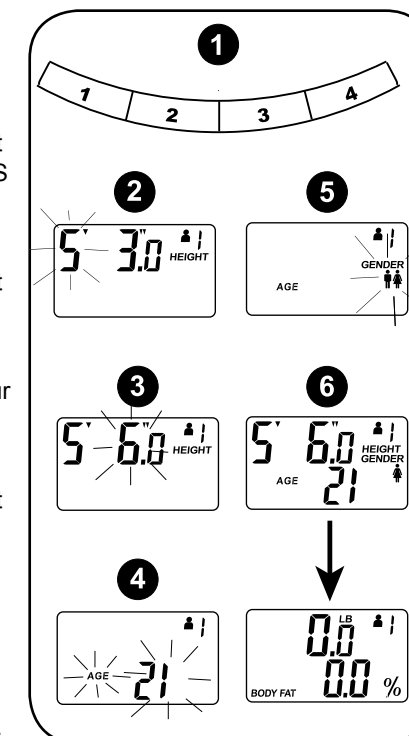Repeat Steps 1-6 to program another user's profile.

NOTE: If a button is not pushed within 20 seconds during programming, the scale will turn off automatically.

If the scale is operating in the Kilograms (KG) mode, the height will be set in centimeters (CM).

Resetting an individual user's history Data

Press and hold the desired User Button and the SET button simultaneously until the HEIGHT icon starts flashing. The user's data has been erased.

5-YEAR LIMITED WARRANTY

PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THE ADDRESSES SHOWN OR TO THE PLACE OF PURCHASE.

Sunbeam Products, Inc. ("Sunbeam") warrants that for a period of five years from the date of purchase, this product will be free from defects in material and workmanship. Sunbeam, at its option, will repair or replace this product or any component of the product found to be defective during the warranty period. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. Sunbeam dealers, service centers, or retail stores selling Sunbeam products do not have the right to alter, modify or in any way change the terms and conditions of this warranty. This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than Sunbeam or an authorized Sunbeam service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

What are the limits on Sunbeam's Liability?

Sunbeam shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition.

Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty.

Sunbeam disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise. Sunbeam shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party.

Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

How to Obtain Warranty Service

If you have any question regarding this warranty or would like to obtain warranty service, please call 1-800-672-5625, or visit us at www.healthometer.com and a convenient service center address will be provided to you.

In the U.S.A, this warranty is offered by Sunbeam Products, Inc. located in Boca Raton, Florida 33431.

PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THE ADDRESSES SHOWN OR TO THE PLACE OF PURCHASE.

©2005 Sunbeam Products, Inc. All Rights Reserved. Health o meter® is a registered trademark of Sunbeam Products, Inc., Boca Raton, FL 33431

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En E.E.U.U., esta garantía es ofrecida por Sunbeam Products, Inc. situada en Boca Raton, Florida 33431.