
# Health o meter®

## WEIGHT & BODY FAT MONITORING SCALE

### User's Manual

Model BFM955

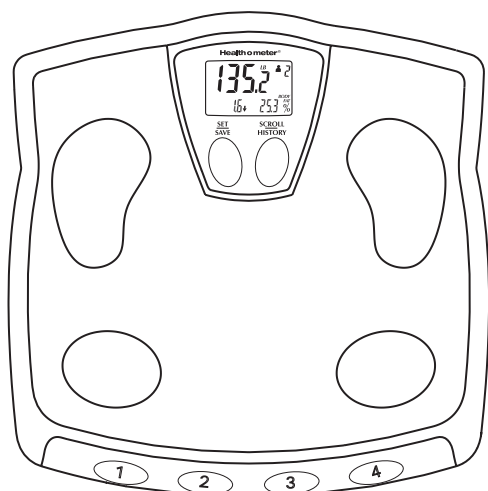
#### WELCOME!

Congratulations on your purchase of this Health o meter® wellness monitoring product. It's so much more than just a scale. It helps monitor your weight and body fat composition for up to four different users. It also tracks weight changes and compares them to your goal. It's easy to use. Once your profile is set, your scale is ready to use. Also, you can save up to five records for each user. It's that simple!


#### B. MONITORING YOUR WEIGHT AND BODY FAT OVER TIME WITHOUT SETTING YOUR GOAL

- Press one of the "USER" buttons. Display will show TIME and DATE; then it will calibrate ("--") and will show "0.0" with either "LB" or "KG" as well as "0.0" body fat.
- Step on the scale and position your feet in a manner that will distribute your weight evenly over the entire scale platform. Stand still while reading the display.
- Step off the scale and press the "SET/SAVE" button to save your reading.

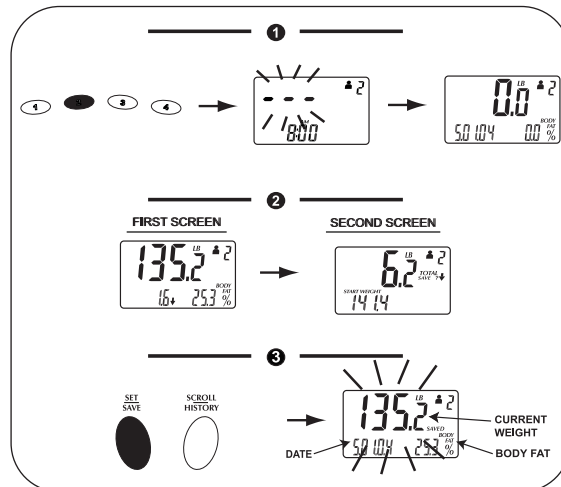
#### C. WEIGHING WITHOUT SELECTING A USER BUTTON

- Tap on the scale platform. The scale will calibrate ("--"), and then will display "0.0" with either "LB" or "KG".
- Step on the scale and position your feet in a manner that will distribute your weight evenly over the entire scale platform. The LCD display will flash "--" several times for approximately 3 seconds before displaying your weight. Stand still while reading the display.

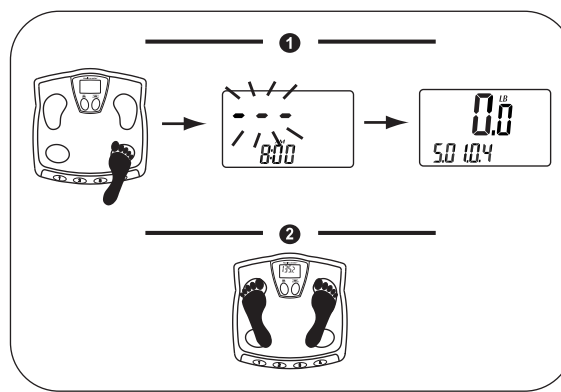**CAUTION** Weight exceeding the capacity (350 lb/159 kg) may damage your scale. This product is designed for personal use only. Not for commercial or industrial use. This unit is not waterproof: avoid contact with excessive moisture.

#### KEYS AND ICONS

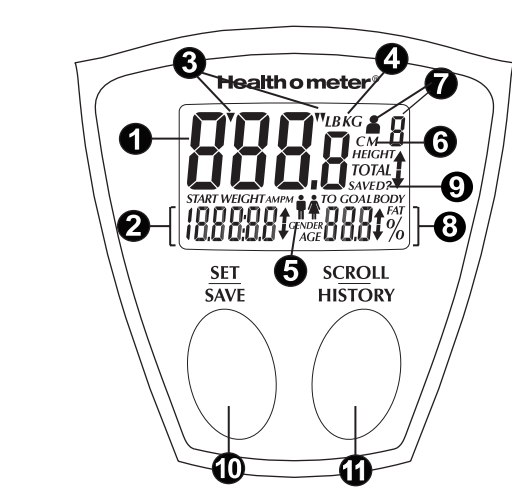- Weight/Height Readout
- Time/Date, start weight and difference in weight compared to previously saved weight
- Height displayed in feet and inches
- Unit of Measure: Pounds (LB) or Kilograms (KG)
- Gender Symbols
- Height displayed in centimeters
- User number and memory indicator
- Body Fat %, Age Readout, Goal, Difference in weight compared to goal and change from previously saved weight
- Information ready to be saved
- SET/SAVE button
- SCROLL/HISTORY button
- User buttons

#### IMPORTANT!

Please read all instructions to familiarize yourself with features and operations before programming the scale.


#### CAUTION!

- The use of this scale is not recommended for persons with an electric implant (heart pacemaker).
- The scale will not calculate body fat unless you have bare feet.
- For an accurate body fat estimate, do not bend your knees, and keep your legs/thighs apart and straight when standing on the scale.
- The body fat function is not recommended for use by pregnant women, persons with fever, swollen legs or other edemas, as well as a person who is over hydrated or dehydrated.
- Body fat percentages may be slightly higher or lower for children under 16, persons with diabetes and other health conditions.
- Women naturally store approximately 5% more body fat than men.

#### REVIEWING SAVED RECORDS (HISTORY)

NOTE: Two screens will be displayed for each record.

- Press the desired "USER" button.
- Press the "SCROLL/HISTORY" button to scroll down the list of saved records.
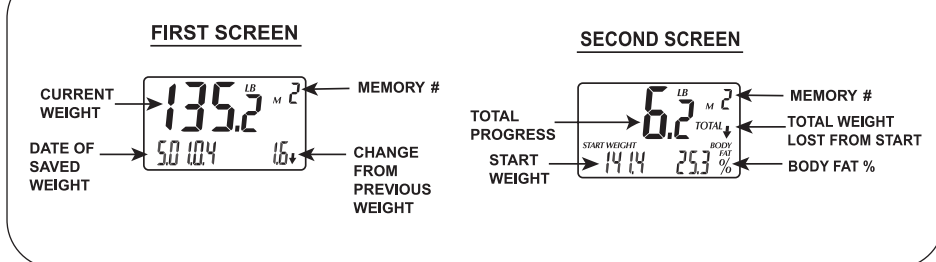
#### RESETTING AN INDIVIDUAL USER'S HISTORY DATA

This process will erase the records for the selected user.

NOTE: You will be prompted to reset your profile and goal.

- Press the "SET/SAVE" button and a "USER" button at the same time for 3 seconds. All your history will be erased. The scale will prompt you to reset your profile. For further instructions, please see page 3.

#### COMPLETELY RESETTING THE SCALE

All data will be lost. Erased data cannot be recovered.

Make sure the LCD display is turned off.

- Turn the scale over.
- Using a pin, press the button inside the small hole located near the top for 3 seconds. The scale's memory has been reset.

#### ABOUT BODY FAT

The Health o meter® Body Fat Monitoring Scale can help you in achieving your physical best by providing two critical measurement estimates: BODY FAT COMPOSITION and BODY WEIGHT. In addition it helps you track your progress to a goal by allowing you to set your weight goal and calculating the weight differences between your starting point, your previously saved weight and your set weight goal.


Traditional bathroom scales, indicate your overall weight—both lean mass and body fat together. The Health o meter® Body Fat Monitoring Scale lets you track your weight and body fat separately. See enclosed brochure for more information on Body Fat.

#### FIRST TIME SETUP


##### Back of the Scale

- Select pounds (LB) or kilograms (KG)
- Place the scale on a flat surface. A hard surface is best for maximum accuracy.

#### ADJUSTING TIME AND DATE

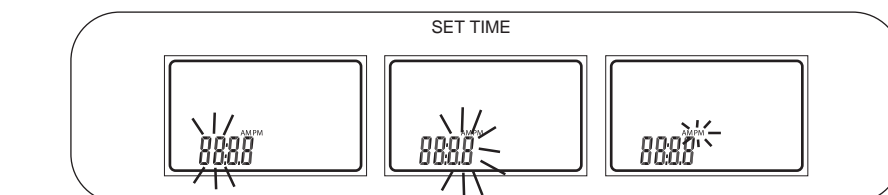- Hold down the "SET/SAVE" button until the HOUR starts blinking.
- Press or hold down the "SCROLL/HISTORY" button to choose the desired HOUR.
- Press the "SET/SAVE" button. MINUTES will start blinking.
- Press or hold down the "SCROLL/HISTORY" button to choose the desired MINUTES.
- Press the "SET/SAVE" button. "AM" or "PM" will start blinking.
- Press the "SCROLL/HISTORY" button to choose "AM" or "PM".
- Press the "SET/SAVE" button. Now your clock is set.

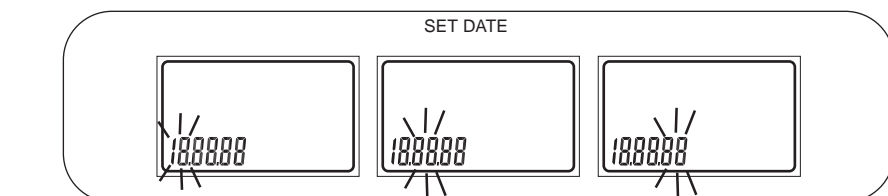
#### REPLACING THE BATTERIES

This unit includes four (4) AA batteries. When "Lo" appears on the display, replace all batteries.

**CAUTION** Removing the batteries will reset the Date and Time but the recorded data will not be affected.

##### Battery Replacement

- Locate the battery cover on the back of the scale. Push the locking tab and pull the battery cover up.
- Use the tip of a pen to loosen each battery. Hold the scale at an angle to allow the batteries to slide out.
- Insert four new batteries, pressing them firmly into their slot.
- Replace the battery cover.
- Set the Time/Date (see ADJUSTING TIME AND DATE).
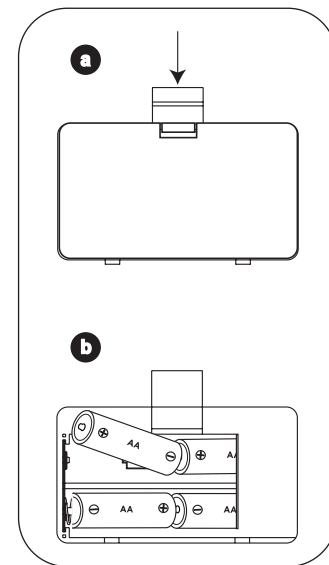
#### TROUBLESHOOTING

If "E" appears on the upper row of the display:

- You may have stepped on the scale before it displayed "0.0".
- The scale cannot detect a stable weight. Stand still while weighing.
- The load has exceeded the scale's capacity (350 lb/159 kg).

If "E" appears on the lower row of the display:

- The scale cannot measure body fat.

NOTE: Display will turn off automatically in 20 seconds if no key is pressed.

**IT IS NORMAL FOR YOUR WEIGHT TO VARY DURING THE DAY AND FROM ONE DAY TO ANOTHER.**

#### CARING FOR YOUR SCALE

Periodically wipe the scale surface down with a DRY cloth to remove dust.

- Press or hold down the "SCROLL/HISTORY" button to choose the desired MONTH.
- Press the "SET/SAVE" button.
- Press or hold down the "SCROLL/HISTORY" button to choose the desired DAY.
- Press the "SET/SAVE" button.
- Press or hold down the "SCROLL/HISTORY" button to choose the desired YEAR.
- Press the "SET/SAVE" button. You have setup time and date.

#### PROGRAMMING / RESETTING YOUR PROFILE

- Press and hold the desired user button until the HEIGHT icon starts flashing.
- Press the "SCROLL/HISTORY" button to choose the desired FEET (or CENTIMETERS in KG mode). Then press the "SET/SAVE" button.
- Press the "SCROLL/HISTORY" button to choose the desired INCHES. Then press the "SET/SAVE" button.

NOTE: Step 3 is not required for KG mode.

Now you are ready to set your age and gender:

- Press or hold down the "SCROLL/HISTORY" button to choose the desired AGE. Then press the "SET/SAVE" button.
- Press the "SCROLL/HISTORY" button to choose the desired gender (MALE or FEMALE). Then press the "SET/SAVE" button.

Display should show HEIGHT, AGE and GENDER. Then "YES" will start blinking.

Now you are ready to set your weight goal:


- Press the "SCROLL/HISTORY" button to select "YES" or "NO".

NOTE: You can set a goal (by selecting "YES") or just track your weight (by selecting "NO").

- Press the "SET/SAVE" button to store your selection.


- Now enter your desired goal weight. Scroll to choose your desired weight and press the "SET/SAVE" to store it.

NOTE: Each digit is set separately, until you have programmed your goal weight.

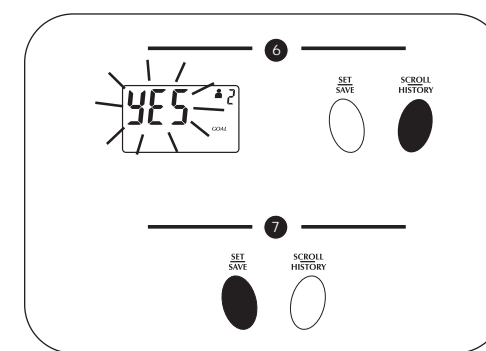
#### BODY FAT COMPOSITION RANGES

##### Men's Body Fat Composition Ranges

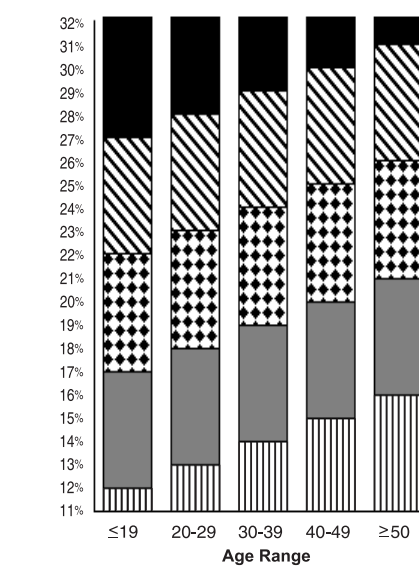
##### Women's Body Fat Composition Ranges

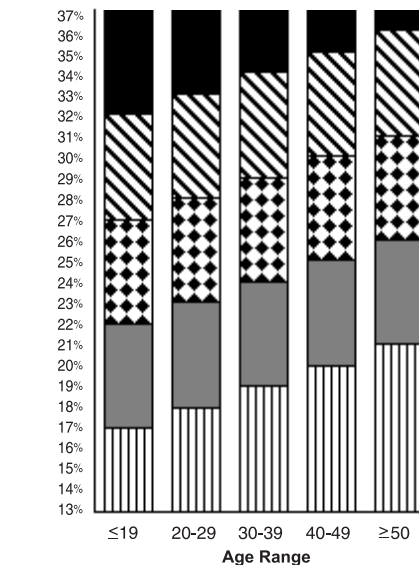
Extracted from Principals & Labs for Fitness and Wellness, 1st Edition ©1999. To determine the percentage of body fat that is appropriate for your body, consult your physician.

- The display will flash "SAVED" twice, followed by "--". Then the scale will calibrate to "0.0" with either "LB" or "KG".
- Your profile and goal (if desired) has been saved and the scale is ready for use.
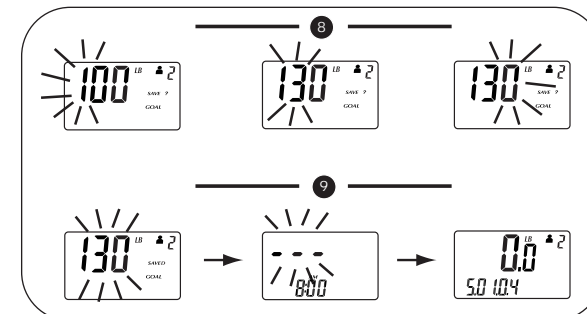
#### USING YOUR BODY FAT MONITORING SCALE

##### A. MONITORING YOUR WEIGHT, BODY FAT PERCENTAGE AND COMPARING IT TO YOUR GOAL.

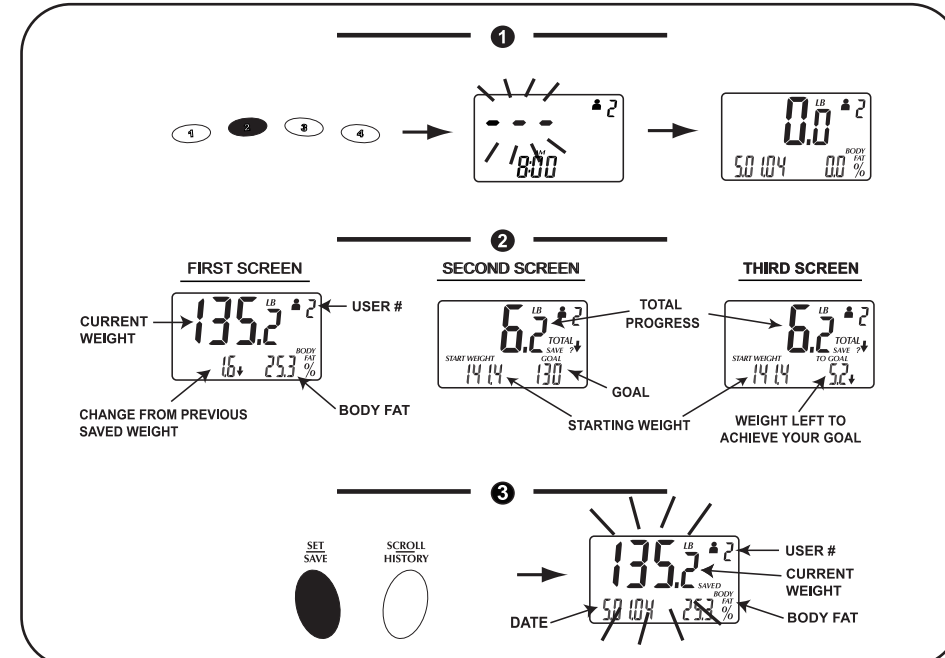- Press one of the "USER" buttons. Display will show TIME, DATE and user profile; then it will calibrate ("--") and will show "0.0" with either "LB" or "KG" as well as "0.0" body fat.
- Step on the scale and position your feet in a manner that will distribute your weight evenly over the entire scale platform. Stand still while reading the display. Three screens will be displayed.
- Step off the scale and press the "SET/SAVE" button to save your reading.

#### 5 YEAR LIMITED WARRANTY

**PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THE ADDRESSES SHOWN OR TO THE PLACE OF PURCHASE.**

Sunbeam Products, Inc. ("Sunbeam") warrants that for a period of five years from the date of purchase, this product will be free from defects in material and workmanship. Sunbeam, at its option, will repair or replace this product or any component of the product found to be defective during the warranty period. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. This warranty is valid for the original retail purchaser from the date of initial purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. Sunbeam dealers, service centers, or retail stores selling Sunbeam products do not have the right to alter, modify or in any way change the terms and conditions of this warranty. This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than Sunbeam or an authorized Sunbeam service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

##### What are the limits on Sunbeam's Liability?

Sunbeam shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition. Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty.

Sunbeam disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise. Sunbeam shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party.

Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

##### How to Obtain Warranty Service

If you have any question regarding this warranty or would like to obtain warranty service, please call 1-800-672-5625, or visit us at [www.healthometer.com](http://www.healthometer.com) and a convenient service center address will be provided to you.

In the U.S.A, this warranty is offered by Sunbeam Products, Inc. located in Boca Raton, Florida 33431.

**PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THE ADDRESSES SHOWN OR TO THE PLACE OF PURCHASE.**

© 2004 Sunbeam Products, Inc. All rights reserved. Health o meter® is a trademark of Sunbeam Products, Inc. Boca Raton, FL 33431. Printed in China.

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En E.E.U.U., esta garantía es ofrecida por Sunbeam Products, Inc. situada en Boca Raton, Florida 33431.