
# Health o meter™

## DIGITAL BODY FAT CALIPER

User's Manual  
Model HDSC012-01

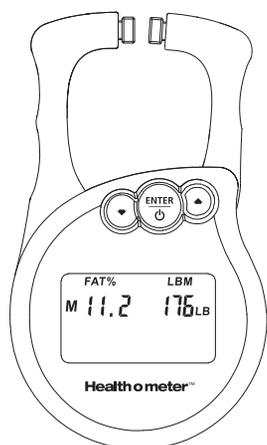www.healthometer.com

### WELCOME!

Congratulations on your purchase of the Health o meter™ digital body fat caliper. This product will estimate your body fat percentage based off of three skin fold sites.


### PLEASE READ CAREFULLY AND SAVE

For information on service or warranty, contact Sunbeam Products, Inc. at 1-800-672-5625 or go to www.healthometer.com.

## ABOUT YOUR HEALTH O METER™ BODY FAT CALIPER

The Health o meter™ Digital Body Fat Caliper estimates your body fat percentage based off of three skin fold sites.

### Keys and Icons

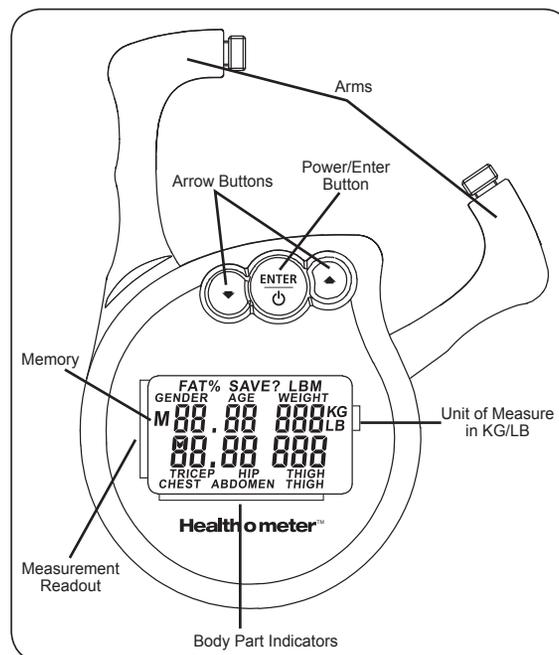**CAUTION** This product is designed for personal use only. Not for commercial or industrial use. This unit is not waterproof; avoid contact with excessive moisture.

## USER PROFILE SET UP

1. Press and hold the **Power/Enter** button to turn the unit on.
2. The metric unit will begin blinking.
3. Press the **up** (▲) or **down** (▼) arrows to select the unit of measurement (lb or kg) and press the **Power/Enter** button to confirm.
4. The weight will begin blinking.
5. Press the **up** (▲) or **down** (▼) arrows to display your weight and press the **Power/Enter** button to confirm.
6. The age number will begin blinking.
7. Press the **up** (▲) or **down** (▼) arrows to display your age and press the **Power/Enter** button to confirm.
8. The gender label will begin blinking.
9. Press the **up** (▲) or **down** (▼) arrows to select male (M) or female (F) and press the **Power/Enter** button to confirm.

The user profile is now set up and will be displayed the next time the unit is turned on. You will also have the ability to change the user profile.

### Turning the unit on and off:

To turn the unit on or off, press and hold the **Power/Enter** button.

### Prior Measurements:

This product will display the prior confirmed measurements. New measurements that are confirmed will overwrite the prior measurement.

### Updating profile or Changing Units of Measurement:


User profile information is automatically saved after each use. Changes in weight, or new individuals wanting to use the product, will require that the profile information be updated accordingly for accurate calculations.

To update the user profile or unit of measurement, simply turn the unit on and update the user information and/or metric as needed. Refer to the User Profile Set Up section.

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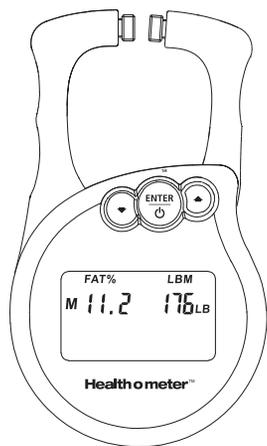## Measuring sites:


**NOTE:** In order to calculate your body fat %, you will need to have your user profile set up prior to measuring the skin fold sites. See User Profile Set Up section.

1. After your user profile has been set up, you can begin using the product to measure the skin fold sites.
2. The product will prompt which skin fold site to be measured by blinking on and off.
3. To get a good skin fold pinch, use your thumb and forefinger to pinch the site. The goal is to pinch deep enough to capture the fat.
4. While holding the pinch, use the skin fold caliper to measure the site. The measurement should be made close to the pinch and about half way up from the base of skin fold.
5. Squeeze the arms of the caliper together to measure the pinch. The unit will make an audible beep when a measurement has been taken.
6. Press the **Power/Enter** button to confirm.
7. The product will cycle to the next skin fold site.
8. Repeat the same steps for each skin fold site.
9. After you have measured and confirmed the three skin fold sites, the unit will automatically calculate your body fat % and lean body mass (body weight minus body fat).
10. The unit will then prompt you to save the information. Press the **Power/Enter** button to confirm. Otherwise, power the unit off by pressing and holding the **Power/Enter** button.


In order for the body fat percentage to be calculated, you will need to measure all three skin fold sites displayed on the bottom of screen.

See images and content below for more information on the specific skin fold sites. Also, pay attention to the angles of the different pinch sites.

### MEN: 3 Skin Fold Sites
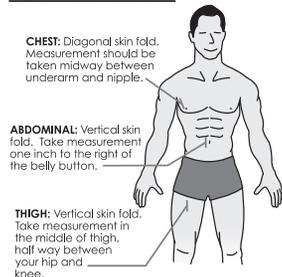
### WOMEN: 3 Skin Fold Sites

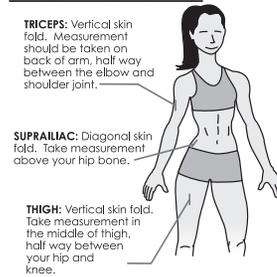Using consistent pinch points is important. Results may vary depending on pinch points.

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
## REPLACING THE BATTERY

1. Locate the battery cover on the back of the tape measure.
2. Use a Phillips head screw driver to remove the screw from the battery cover.
3. Remove the used battery and insert a new Lithium battery (CR2032-3V).
4. Replace the battery cover and secure the screw with Phillips head screw driver.

**NOTE:** The Lithium coin battery(s) used in this product contains perchlorate material. Special handling may apply in California. Go to: [www.dtsc.ca.gov/hazardouswaste/perchlorate](http://www.dtsc.ca.gov/hazardouswaste/perchlorate) for more information.
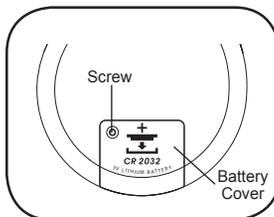
## CARING FOR YOUR BODY FAT CALIPER

Periodically wipe the surface of the tape measure with a DRY cloth to remove dust.

## POTENTIAL FOR RADIO/TELEVISION INTERFERENCE


This product has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

The product generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that the interference will not occur in a particular installation. If the product does cause harmful interference to radio or television reception, which can be determined by turning the product on or off, the user is encouraged to try to correct the interference by one or more of the following measures: a) Reorient or relocate the receiving antenna; b) Increase the separation between the product and the receiver; c) Connect the product into an outlet on a circuit different from that to which the receiver is connected; d) Consult the dealer or an experienced radio/TV technician for help. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

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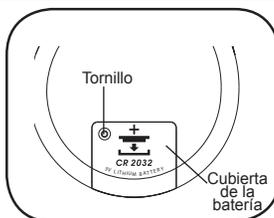## 5 YEAR LIMITED WARRANTY

**PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE.**

Sunbeam Products, Inc. doing business as Jarden Consumer Solutions or if in Canada, Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions (collectively "JCS") warrants that for a period of five years from the date of purchase, this product will be free from defects in material and workmanship. JCS, at its option, will repair or replace this product or any component of the product found to be defective during the warranty period. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. Do NOT attempt to repair or adjust any electrical or mechanical functions on this product. Doing so will void this warranty. This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. JCS dealers, service centers, or retail stores selling JCS products do not have the right to alter, modify or in any way change the terms and conditions of this warranty. This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than JCS or an authorized JCS service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

### What are the limits on JCS's Liability?

JCS shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition. Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty. JCS disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise. JCS shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party. Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you. This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

### How to Obtain Warranty Service

**In the U.S.A.:** If you have any questions regarding this warranty or would like to obtain warranty service, please call 1 800-672-5625 and a convenient service center address will be provided to you.

**In Canada:** If you have any questions regarding this warranty or would like to obtain warranty service, please call 1 800-667-8623 and a convenient service center address will be provided to you.

**In the U.S.A.,** this warranty is offered by Sunbeam Products, Inc. doing business as Jarden Consumer Solutions located in Boca Raton, Florida 33431. In Canada, this warranty is offered by Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions, located at 20 B Hereford Street, Brampton, Ontario L6Y 0M1. If you have any other problem or claim in connection with this product, please write our Consumer Service Department.

**PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE.**

© 2012 Sunbeam Products, Inc. doing business as Jarden Consumer Solutions. All rights reserved. Distributed by Sunbeam Products, Inc. doing business as Jarden Consumer Solutions, Boca Raton, Florida 33431. Printed in China.

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