Massaging Renue™ XL Tension Relieving Heat Therapy Pad

Model 882

Congratulations! You have purchased Massaging Renue™ XL Tension Relieving Heat Therapy Pad. Sit back, relax and start enjoying the soothing benefits of the Massaging Renue™ XL Tension Relieving Heat Therapy Pad.

# IMPORTANT SAFETY INSTRUCTIONS DO NOT DESTROY

**DANGER** – TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK AND FIRE, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

1. READ ALL INSTRUCTIONS.
2. DO NOT USE WHILE SLEEPING.
3. DO NOT USE ON AN INFANT OR ON AN ANIMAL.
4. THIS PAD IS NOT TO BE USED ON OR BY AN INVALID, SLEEPING, HELPLESS OR UNCONSCIOUS PERSON, A PERSON WITH DIABETES, A PERSON WITH POOR BLOOD CIRCULATION OR INSENSITIVITY TO HEAT.
5. THIS APPLIANCE IS NOT TO BE USED BY CHILDREN. CLOSE SUPERVISION IS NECESSARY WHEN ANY APPLIANCE IS USED NEAR CHILDREN. CHILDREN SHOULD NOT PLAY WITH THE APPLIANCE.
6. DO NOT USE ON AREAS OF INSENSITIVE SKIN.
7. BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING, CHECK SKIN UNDER PAD FREQUENTLY.
8. IF PROVIDED WITH A REMOVABLE COVER, NEVER USE PAD WITHOUT COVER IN PLACE.
9. DO NOT USE IN AN OXYGEN ATMOSPHERE.
10. DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.
11. DO NOT SIT ON, AGAINST, OR CRUSH PAD – AVOID SHARP FOLDS. PLACE PAD ON TOP OF AND NOT UNDER THE PART OF BODY NEEDING HEAT.
12. UNPLUG WHEN NOT IN USE. NEVER LEAVE APPLIANCE UNATTENDED, ESPECIALLY IF CHILDREN ARE PRESENT.
13. NEVER PULL THIS PAD BY THE SUPPLY CORD. DO NOT USE THE CORD AS A HANDLE. DO NOT PINCH CORD.
14. DO NOT TUCK IN, TRAP, FOLD, CROSS, PINCH OR WRAP CONTROL OR CORD WITH HEATING PAD DURING USE.
15. CAREFULLY EXAMINE PAD, COVER, CORD AND CONTROLLER BEFORE EACH USE. DISCARD THE PAD IF ANY SIGN OF DETERIORATION, WEAR OR DAMAGE. DO NOT TAMPER WITH THIS PAD IN ANY WAY. THERE ARE NO USER SERVICEABLE PARTS. IF PAD DOES NOT FUNCTION SATISFACTORILY, SEE WARRANTY FOR

## CUSTOMER CONTACT INFORMATION.

16. LOOP CORD LOOSELY WHEN STORING, TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS. ALLOW TO COOL BEFORE STORING AND DO NOT CREASE BY PLACING ITEMS ON TOP OF IT. STORE IN THE ORIGINAL PACKAGING, IN A COOL, DARK, DRY PLACE.
17. IF PROVIDED WITH A DETACHABLE CONTROLLER, USE ONLY WITH THE PROVIDED CONTROLLER.
18. USE THIS PAD ONLY AT RATED VOLTAGE, AS MARKED ON THE PRODUCT. DO NOT USE WITH GENERATORS, POWER CONVERTERS, OR INVERTERS.
19. DO NOT USE THIS PAD WITH LINIMENT, SALVE OR OINTMENT PREPARATIONS. SKIN BURNS COULD RESULT.
20. IF PAD IS PROVIDED WITH MAGNETIC CLOSURES: DO NOT USE IF YOU HAVE A PACEMAKER, ICD (IMPLANTABLE CARDIOVERTER DEFIBRILLATOR), INTERNAL INSULIN PUMP OR OTHER MEDICAL DEVICES SENSITIVE TO MAGNETS.
21. THIS APPLIANCE IS NOT INTENDED FOR MEDICAL USE IN HOSPITALS.
22. MODELS WITH A MESSAGE FUNCTION: DO NOT MESSAGE ABDOMEN WHEN PAIN IS PRESENT, OR AREAS THAT ARE SWOLLEN, INFLAMED, COVERED WITH SKIN ERUPTIONS OR UNEXPLAINED PAIN. CONSULT A PHYSICIAN BEFORE USING ANY EQUIPMENT OF THIS TYPE.
23. SOME MODELS MAY CONTAIN FILLING MATERIAL (PELLETS OR BEADS), MADE FROM EITHER ABS PLASTIC (NON-MICROWAVABLE MODELS) OR SILICON DIOXIDE (MICROWAVABLE MODELS). IF FILLING MATERIAL IS RELEASED DUE TO PRODUCT DAMAGE, DISCONTINUE USE OF PRODUCT. DO NOT INGEST FILLING MATERIAL. IF FILLING MATERIAL IS INGESTED, SEEK MEDICAL ATTENTION.
24. SAVE THESE INSTRUCTIONS.

# HOUSEHOLD USE ONLY

## NORTH AMERICAN MODELS WITH POLARIZED PLUGS:

**POLARIZED PLUG**  
THIS APPLIANCE HAS A POLARIZED PLUG, (ONE BLADE IS WIDER THAN THE OTHER). AS A SAFETY FEATURE TO REDUCE THE RISK OF ELECTRICAL SHOCK, THIS PLUG IS INTENDED TO FIT IN A POLARIZED OUTLET ONLY ONE WAY. IF THE PLUG DOES NOT FIT FULLY IN THE OUTLET, REVERSE THE PLUG. IF IT STILL DOES NOT FIT, CONTACT A QUALIFIED ELECTRICIAN. DO NOT ATTEMPT TO DEFEAT THIS SAFETY FEATURE OR MODIFY THE PLUG IN ANY WAY. IF THE PLUG FITS LOOSELY INTO THE AC OUTLET OR IF THE AC OUTLET FEELS WARM DO NOT USE THAT OUTLET.

Massaging Renue™ XL Tension Relieving Heat Therapy Pad is specifically designed to fit the neck, shoulders & upper back area, where users experience the most physical signs of tension and stress. Drape the Massaging Renue™ XL Heat Therapy Pad comfortably on your shoulders making sure that the neck flap is properly in place. Weighted edges and adjustable snap enclosures help to secure the Massaging Renue™ XL Heat Therapy Pad in place.

## HOW TO USE YOUR MASSAGING RENUETM XL TENSION RELIEVING HEAT THERAPY PAD

The Massaging Renue™ XL Tension Relieving Heat Therapy Pad has 2 heat settings and 4 massage settings which can be used simultaneously or individually for customized therapy (Figure 1).

1. Plug into 120V AC electrical outlet.
2. To turn on select either Low/ High Heat for heat therapy, or Low/ High Massage for massage therapy.  
**NOTE:** The 4 Massage settings consists of the following speeds which can be used along with the heat therapy options: (1) Low vibration - consistent, (2) Low vibration - pulsing, (3) High vibration - consistent, (4) High vibration - pulsing.
3. Select desired heat and/ or massage level. **To turn OFF, you must scroll through the remaining heat and/ or massage settings by pressing the buttons sequentially.**
4. Allow the Massaging Renue™ XL Heat Therapy Pad to preheat for several minutes. The Massaging Renue™ XL Heat Therapy Pad has been designed to meet the maximum heat output according to UL Standards.

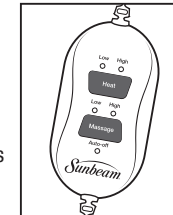
Figure 1

**NOTE:** A red blinking light indicates an error has occurred. Unplug from outlet, verify physical condition of the cord and function of AC outlet. If blinking light persists, discontinue use and contact Consumer Services.

## AUTO-OFF FEATURE

Your Massaging Renue™ XL Heat Therapy Pad controller includes an Auto-Off feature. This feature is designed to conserve energy by automatically turning off after 2 hours of continuous use when using the Heat with Massaging option, and a 30-minute Auto-Off feature when using only the Massage option. The red Auto-Off light will blink, indicating that this feature has activated. **NOTE:** After a power outage, the Auto-Off feature automatically resets the timer.

## WASHING INSTRUCTIONS

**CAUTION: DO NOT DRY CLEAN! DO NOT USE BLEACH! DO NOT USE WRINGER! DO NOT WET INNER PAD CONNECTED WITH CORD! ONLY WASH REMOVABLE OUTER COVER!**

1. Disconnect the Massaging Renue™ XL Heat Therapy Pad from outlet.
2. Remove the outer cover from the Massaging Renue™ XL Heat Therapy Pad.
3. Machine wash the cover ONLY in cold water on gentle cycle, or hand wash.
4. Tumble dry on low heat or line dry.
5. Do not wash heating pad. Wipe heating pad with a damp cloth only, as needed. Do not dampen the cord and connector. Hang to dry.

## FCC CERTIFIED

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

**NOTE:** This product has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This product generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
  - Increase the separation between the product and receiver.
  - Connect the product into an outlet on a circuit different from that to which the receiver is connected.
  - Consult the dealer or an experienced radio/TV technician for help.
- Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

## ICES CERTIFIED

This product complies with CAN ICES-3 (B)/NMB-3(B).

## 2-YEAR LIMITED WARRANTY

Sunbeam Products, Inc. doing business as Jarden Consumer Solutions or if in Canada, Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions (collectively "JCS") warrants that for a period of two years from the date of purchase, this product will be free from defects in material and workmanship. JCS, at its option, will repair or replace this product or any component of the product found to be defective during the warranty period. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. Do NOT attempt to repair or adjust any electrical or mechanical functions on this product. Doing so will void this warranty. This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. JCS dealers, service centers, or retail stores selling JCS products do not have the right to alter, modify or any way change the terms and conditions of this warranty. This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than JCS or an authorized JCS service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

### What are the limits on JCS's Liability?

JCS shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition.

Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty.

JCS disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise.

JCS shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party.

Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

### How to Obtain Warranty Service

**In the U.S.A.** - If you have any question regarding this warranty or would like to obtain warranty service, please call 1-800-435-1250 and a convenient service center address will be provided to you.

**In Canada** - If you have any question regarding this warranty or would like to obtain warranty service, please call 1-800-667-8623 and a convenient service center address will be provided to you.

**In the U.S.A.**, this warranty is offered by Sunbeam Products, Inc. doing business as Jarden Consumer Solutions located in Boca Raton, Florida 33431. **In Canada**, this warranty is offered by Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions, located at 20 B Hereford Street, Brampton, Ontario L6Y 0M1. If you have any other problem or claim in connection with this product, please write our Consumer Service Department.

**PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE.**

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©2017 Sunbeam Products, Inc. operando bajo el nombre de Jarden Consumer Solutions. Distribuido por Sunbeam Products, Inc. operando bajo el nombre de Jarden Consumer Solutions, Boca Raton, Florida 33431. Sunbeam Products, Inc. es una filial de Newell Brands, Inc. (NYSE: NWL).