XpressHeat® Wrapping Pad

For product questions: Sunbeam Consumer Service

USA: 1.800.435.1250

www.sunbeam.com

Models: 802 & 803

IMPORTANT SAFETY INSTRUCTIONS DO NOT DESTROY

DANGER: TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK AND FIRE, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

1. READ ALL INSTRUCTIONS.
2. DO NOT USE WHILE SLEEPING.
3. DO NOT USE ON AN INFANT OR ON AN ANIMAL.
4. THIS PAD IS NOT TO BE USED ON OR BY AN INVALID, SLEEPING, HELPLESS OR UNCONSCIOUS PERSON, A PERSON WITH DIABETES, A PERSON WITH POOR BLOOD CIRCULATION OR INSENSITIVITY TO HEAT, OR A PERSON INCAPACITATED BY MEDICATION.
5. THIS APPLIANCE IS NOT TO BE USED BY CHILDREN. CLOSE SUPERVISION IS NECESSARY WHEN ANY APPLIANCE IS USED NEAR CHILDREN. CHILDREN SHOULD NOT PLAY WITH THE APPLIANCE.
6. DO NOT USE ON AREAS OF INSENSITIVE SKIN.
7. BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING. CHECK SKIN UNDER PAD FREQUENTLY.
8. DO NOT USE IN AN OXYGEN ATMOSPHERE.
9. DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.
10. DO NOT SIT ON, AGAINST, OR CRUSH PAD – AVOID SHARP FOLDS. PLACE PAD ON TOP OF AND NOT UNDER THE PART OF BODY NEEDING HEAT.
11. UNPLUG WHEN NOT IN USE. NEVER LEAVE APPLIANCE UNATTENDED, ESPECIALLY IF CHILDREN ARE PRESENT.
12. NEVER PULL THIS PAD BY THE SUPPLY CORD. DO NOT USE THE CORD AS A HANDLE. DO NOT PINCH CORD.
13. DO NOT TUCK IN, TRAP, FOLD, CROSS, PINCH OR WRAP CONTROL OR CORD WITH HEATING PAD DURING USE.
14. DO NOT PLACE ANYTHING OVER THE HEATING PAD

WHILE IN USE OR PLUGGED IN, SUCH AS A TOWEL, BLANKET, CLOTHING, OR ANY OTHER INSULATING MATERIAL, AND NEVER ALLOW HEATING PAD TO WRAP AROUND ITSELF.

15. CAREFULLY EXAMINE PAD, COVER, CORD AND CONTROLLER BEFORE EACH USE. DISCARD THE PAD IF ANY SIGN OF DETERIORATION, WEAR OR DAMAGE. DO NOT TAMPER WITH THIS PAD IN ANY WAY. THERE ARE NO USER SERVICEABLE PARTS. IF PAD DOES NOT FUNCTION SATISFACTORILY, SEE WARRANTY FOR CUSTOMER CONTACT INFORMATION.
16. LOOP CORD LOOSELY WHEN STORING, TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS. ALLOW TO COOL BEFORE STORING AND DO NOT CREASE BY PLACING ITEMS ON TOP OF IT. STORE IN THE ORIGINAL PACKAGING, IN A COOL, DARK, DRY PLACE.
17. IF PROVIDED WITH A DETACHABLE CONTROLLER, USE ONLY WITH THE PROVIDED CONTROLLER.
18. USE THIS PAD ONLY AT RATED VOLTAGE, AS MARKED ON THE PRODUCT. DO NOT USE WITH GENERATORS, POWER CONVERTERS, OR INVERTERS.
19. DO NOT USE THIS PAD WITH LINIMENT, SALVE OR OINTMENT PREPARATIONS. SKIN BURNS COULD RESULT.
20. THIS APPLIANCE IS NOT INTENDED FOR MEDICAL USE IN HOSPITALS.
21. SAVE THESE INSTRUCTIONS.

HOUSEHOLD USE ONLY

North American models with polarized plugs:



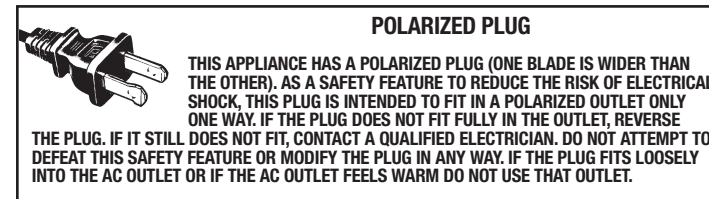
Figure 1
Model 802
Standard Straps

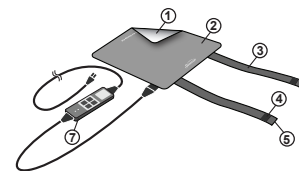1. Interior pad fabric
2. Exterior pad fabric
3. Elastic fastening straps
4. Hook fastener

Model 803
Extendable Straps

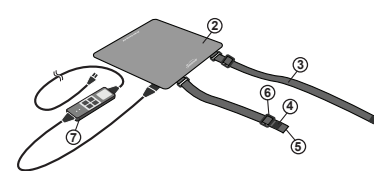5. Strap release tab
6. Plastic strap length adjuster
7. Controller


HOW TO APPLY YOUR HEATING PAD

Note: Always consult with your physician or health professional for recommended treatment of an injury or chronic muscle or joint pain.

The XpressHeat® Wrapping Pad has two elastic fastening straps that allow you to secure the pad in place for hands-free, targeted heat therapy.

1. This heating pad can be worn on several different areas of the body. The elastic fastening straps (Figure 1, number 3) stretch for size adjustability and a comfortable fit. See Figure 2 for ways to wear your heating pad.
2. Drape heating pad around the area requiring heat therapy with interior pad fabric toward body (Figure 1, number 1). Position pad so controller (Figure 1, number 7) is easy to access.
3. Wrap each elastic strap around body and securely fasten directly to the exterior pad fabric (Figure 1, number 2) using the hook fastener (Figure 1, number 4). The hook fastener can also be secured to the strap itself.
CAUTION: Do not allow heated area of pad to overlap onto itself.
4. To remove pad, detach straps by pulling up on the strap release tabs (Figure 1, number 5). Unwrap pad from body.
5. Some wear of the exterior pad fabric may be seen over time; this is normal and will not affect the performance of the heating pad. Use caution when fastening straps to prevent hook fasteners from touching clothing.


Figure 2- APPLICATIONS FOR MODELS 802 & 803

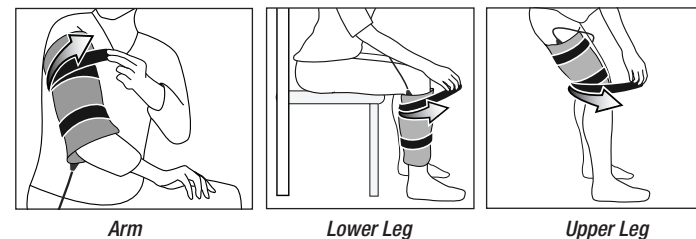ADDITIONAL APPLICATION INSTRUCTIONS FOR EXTENDABLE STRAPS (MODEL 803)

1. Model 803 can be used on larger areas, such as the back and shoulder. See Figure 3. Use the plastic strap length adjuster to extend or retract strap to desired length (Figure 1, number 6).
2. To lengthen or shorten strap length, hold plastic strap length adjuster and pull firmly on the elastic fastening strap.
3. For Back: Extend both straps to appropriate length to fit around waist. Place pad on back and wrap straps around waist securing the hook fasteners to exterior pad fabric (Figure 3).
4. For Shoulder: Adjust straps so that the elastic fastening strap closest to arm is retracted to a shorter length, and strap closest to shoulder is extended to a longer length. Position pad on shoulder with straps hanging down your back. The longer strap should be on the inside of your back and the shorter strap should be on the outside closest to your arm. Wrap shortened strap around arm and secure hook fastener to exterior pad fabric. Then wrap the lengthened strap around back and across the chest to fit pad snugly around shoulder (Figure 3).

Figure 3- ADDITIONAL APPLICATIONS FOR EXTENDABLE STRAPS (MODEL 803)

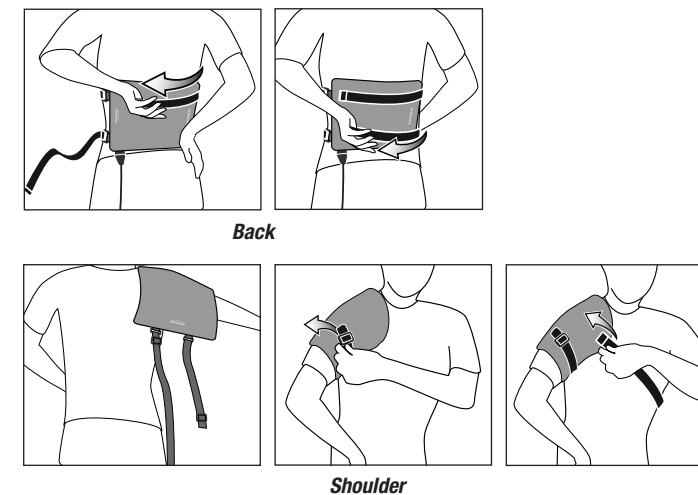HOW TO USE YOUR HEATING PAD

1. Make sure Cord Connector is securely attached to the Pad Connector (Figure 4). Push both parts together for a tight fit.
NOTE: Control must be connected to heating pad before plugging unit into outlet.
2. After control is connected to pad, plug into 120V AC electrical outlet.
3. To turn the pad on, press the Power Button (Figure 5, letter B).
4. Select from 6 heat settings by scrolling to the desired setting with HEAT▲ or HEAT▼ buttons (Figure 5, letter F).
a. Your Sunbeam® XpressHeat® Wrapping Pad features a unique fast-heat technology in which you will feel heat in 30 seconds - for faster pain relief!
The light ring around the controls, called the Heat Meter (Figure 5, letter E), will light up one side at a time, indicating progress as the pad heats up. All four sides of the Heat Meter will light up within 30 seconds of powering up, indicating that the heating pad is warm and ready for use.
5. To turn the pad OFF, press the Power Button one more time. The heating pad and lights of the Heat Meter will turn off.

- CAUTION:**
- DO NOT SIT ON, LEAN AGAINST, CRUSH OR FOLD PAD DURING USE.
 - DO NOT USE WHILE SLEEPING.
 - SEE ADDITIONAL WARNINGS IN THE IMPORTANT SAFETY INSTRUCTIONS SECTION.

AUTO-OFF FEATURE

Your heating pad controller includes a selectable Auto-Off feature. By pressing the "Mode" button (Figure 5, letter C), you can set the timer to "Stay On" continuously or to "2 Hour Auto Off" to shut the pad down. The Auto Off Indicator (Figure 5, letter D) will light up to show selected mode. **NOTE:** After a power outage, the Auto-Off feature automatically resets the timer.

HOW TO USE YOUR HEATING PAD WITH MOIST HEAT THERAPY

In some cases, you may find that moist heat therapy provides more comfort and relief. *When using moist heat therapy, you should protect furniture and clothing from becoming damp.*

1. Dampen the interior pad fabric (Figure 1, number 1) of the heating pad with water by patting the surface with a wet sponge or cloth or use a spray bottle to apply a fine mist. Do not thoroughly soak pad.
2. Do not spray or dampen the pad connector. (Figure 6)

DANGER: DO NOT IMMERSE THE PAD OR WIRING IN WATER OR PLACE UNDER RUNNING WATER. FAILURE TO COMPLY MAY RESULT IN PRODUCT FAILURE AND/OR PERSONAL INJURY.

THE ELECTRICAL PAD CONNECTOR MUST BE THOROUGHLY DRY BEFORE ESTABLISHING ELECTRICAL CONNECTION.

WASHING INSTRUCTIONS

CAUTION: DO NOT DRY CLEAN! DO NOT BLEACH! DO NOT USE WRINGER!

1. Disconnect pad from outlet. DO NOT PULL BY THE CORD.
2. Disconnect controller cord from pad by holding the pad side connector with one hand and pulling firmly on the cord side connector with the other hand.
3. Machine wash pad in cold water on gentle cycle or hand wash. Tumble dry on low heat and remove promptly from dryer or drape to dry.
4. DO NOT plug in pad and pad connector until completely dry.

Figure 4


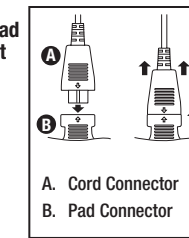
Figure 5

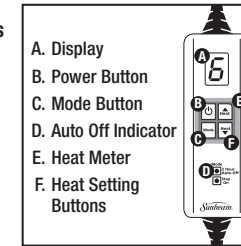FCC CERTIFIED

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This product has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This product generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the product and receiver.
- Connect the product into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Figure 6

5-YEAR LIMITED WARRANTY

Sunbeam Products, Inc. warrants that for a period of five years from the date of purchase, this product will be free from defects in material and workmanship. Sunbeam, at its option, will repair or replace this product or any component of the product found to be defective during the warranty period. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. Do NOT attempt to repair or adjust any electrical or mechanical functions on this product. Doing so will void this warranty. This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. Sunbeam dealers, service centers, or retail stores selling Sunbeam products do not have the right to alter, modify or in any way change the terms and conditions of this warranty. This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than Sunbeam or an authorized Sunbeam service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornados.

What are the limits on Sunbeam's Liability?

Sunbeam shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition. Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty.

Sunbeam disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise. Sunbeam shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party. Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you. This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

How to Obtain Warranty Service

In the U.S.A. - If you have any question regarding this warranty or would like to obtain warranty service, please call 1-800-435-1250 and a convenient service center address will be provided to you.

PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE.

For inquiries regarding recycling and proper disposal of this product, please contact your local waste management facility.

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