For product questions:  
Sunbeam Consumer Service  
USA: 1.800.435.1250  
[www.sunbeam.com](http://www.sunbeam.com)

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## FlexFit™ Heated Wrap Almohadilla Térmica Envoltoria FlexFit™

Instruction Manual  
Manual de Instrucciones

© 2019 Sunbeam Products, Inc. All rights reserved. Distributed by Sunbeam Products, Inc., Boca Raton, Florida 33431. Sunbeam Products, Inc. is a subsidiary of Newell Brands Inc. (NYSE: NWL).

© 2019 Sunbeam Products, Inc. Todos los derechos reservados. Distribuido por Sunbeam Products, Inc., Boca Raton, Florida 33431. Sunbeam Products, Inc. es una subsidiaria de Newell Brands Inc. (NYSE: NWL).

FlexFit 900\_19ESM2 (US)  
Printed in China / Impreso en China

P.N. 196831 Rev A  
GCDS-SUN476667-LD

Model / Modelo: 900

# **IMPORTANT SAFETY INSTRUCTIONS.**

## **DO NOT DESTROY.**

**DANGER:** TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK AND FIRE, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

1. READ ALL INSTRUCTIONS.
2. DO NOT USE WHILE SLEEPING.
3. DO NOT USE ON AN INFANT OR ON AN ANIMAL.
4. THIS PAD IS NOT TO BE USED ON OR BY AN INVALID, SLEEPING, HELPLESS OR UNCONSCIOUS PERSON, A PERSON WITH DIABETES, A PERSON WITH POOR BLOOD CIRCULATION OR INSENSITIVITY TO HEAT, OR A PERSON INCAPACITATED BY MEDICATION.
5. THIS APPLIANCE IS NOT TO BE USED BY CHILDREN. CLOSE SUPERVISION IS NECESSARY WHEN ANY APPLIANCE IS USED NEAR CHILDREN. CHILDREN SHOULD NOT PLAY WITH THE APPLIANCE.
6. DO NOT USE ON AREAS OF INSENSITIVE SKIN.
7. DO NOT OVERLAP THE HEATED AREA LOCATED

- BETWEEN THE RED EDGES. BURNS MAY OCCUR.
8. BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING. CHECK SKIN UNDER PAD FREQUENTLY.
9. DO NOT USE IN AN OXYGEN ATMOSPHERE.
10. DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.
11. DO NOT WRAP PAD TOO TIGHT AROUND THE BODY. PAD SHOULD BE SNUG, BUT SHOULD NOT CUT OFF CIRCULATION.
12. DO NOT WRAP PAD AROUND OR NEAR THE NECK TO AVOID PERSONAL INJURY.
13. DO NOT SIT ON, AGAINST OR CRUSH PAD – AVOID SHARP FOLDS. PLACE PAD ON TOP OF AND NOT UNDER THE PART OF BODY NEEDING HEAT.
14. UNPLUG WHEN NOT IN USE. NEVER LEAVE APPLIANCE UNATTENDED, ESPECIALLY IF CHILDREN ARE PRESENT.
15. NEVER PULL THIS PAD BY THE SUPPLY CORD. DO NOT USE THE CORD AS A HANDLE. DO NOT PINCH CORD.
16. DO NOT TUCK IN, TRAP, FOLD, CROSS, PINCH, OR WRAP CONTROL OR CORD WITH HEATING PAD DURING USE.
17. USE CAUTION WITH PAD ENDS TO PREVENT HOOK FASTENERS FROM DAMAGING CLOTHING OR OTHER

**PROPERTY.**

18. **DO NOT PLACE ANYTHING OVER THE HEATING PAD WHILE IN USE OR PLUGGED IN, SUCH AS A TOWEL, BLANKET, CLOTHING, OR ANY OTHER INSULATING MATERIAL, AND NEVER ALLOW HEATING PAD TO WRAP AROUND ITSELF.**
19. **CAREFULLY EXAMINE PAD, COVER, CORD, AND CONTROLLER BEFORE EACH USE. DISCARD THE PAD IF ANY SIGN OF DETERIORATION, WEAR OR DAMAGE. DO NOT TAMPER WITH THIS PAD IN ANY WAY. THERE ARE NO USER SERVICABLE PARTS. IF PAD DOES NOT FUNCTION SATISFACTORILY, SEE WARRANTY FOR CUSTOMER CONTACT INFORMATION.**
20. **LOOP CORD LOOSELY WHEN STORING, TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS. ALLOW TO COOL BEFORE STORING AND DO NOT CREASE BY PLACING ITEMS ON TOP OF IT. STORE IN THE ORIGINAL PACKAGING, IN COOL, DARK, DRY PLACE.**
21. **IF PROVIDED WITH A DETACHABLE CONTROLLER, USE ONLY WITH THE PROVIDED CONTROLLER.**
22. **USE THIS PAD ONLY AT RATED VOLTAGE, AS MARKED ON THE PRODUCT. DO NOT USE WITH GENERATORS, POWER CONVERTERS, OR**

**INVERTERS.**

23. **DO NOT USE THIS PAD WITH LINIMENT, SALVE OR OINTMENT PREPERATIONS. SKIN BURNS COULD RESULT.**
24. **THIS APPLIANCE IS NOT INTENDED FOR MEDICAL USE IN HOSPITALS.**
25. **SAVE THESE INSTRUCTIONS.**

**North American models with polarized plugs:**

**POLARIZED PLUG**


This appliance has a polarized plug, (one blade is wider than the other). As a safety feature to reduce the risk of electrical shock, this plug is intended to fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature or modify the plug in any way. If the plug fits loosely into the AC outlet or if the AC outlet feels warm do not use that outlet.

**HOUSEHOLD USE ONLY**

## HEATING PAD COMPONENTS

Figure 1

Model: 900

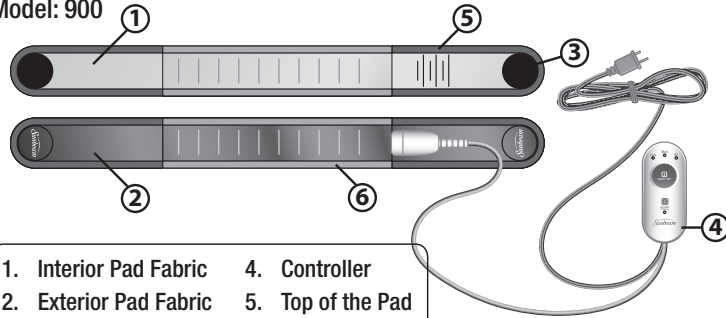
- |                        |                   |
|------------------------|-------------------|
| 1. Interior Pad Fabric | 4. Controller     |
| 2. Exterior Pad Fabric | 5. Top of the Pad |
| 3. Hook Fastener       | 6. Red Edges      |

## WHERE TO USE YOUR HEATING PAD

Heating pads are effective pain relievers to soothe muscle and joint pain. The FlexFit™ Heated Wrap has a long slim design with fastening ends that allow you to wrap the heating pad to multiple areas of the body and secure the pad in place for hands-free, targeted heat therapy.

This heating pad can be worn on several different areas of the body to customize treatment. It is ideal for use on the upper and lower arm, wrist, lower leg, and ankle.


Figure 2
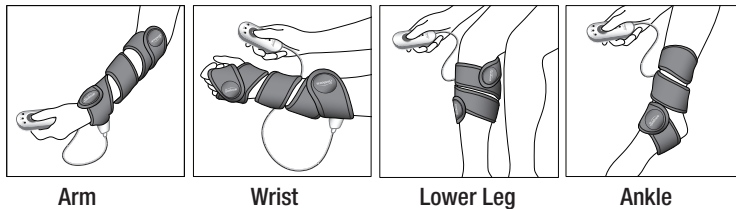
## HOW TO APPLY YOUR HEATING PAD

**Note:** Always consult with your physician or health professional for recommended treatment of an injury or chronic muscle or joint pain.

Figure 3

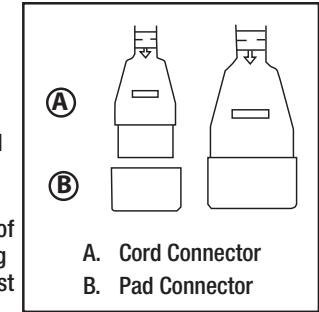
1. Ensure the heating pad is connected to the controller by making sure the cord connector is securely attached to the pad connector before using (Figure 3). Push both parts together for a tight fit.
2. To wrap the heating pad, start at the top of the pad where it indicates 'start wrapping here' (Figure 1, number 5). Secure the first hook fastener (Figure 1, number 3). This is the end where the pad connects to the controller.
3. Wrap the pad around the body and securely fasten the second hook fastener directly to the exterior pad fabric. See figure 2 for the best way to wrap each body part. **NOTE:** Be careful not to overlap the pad. The red edges of the pad should never touch the pad (Figure 1, number 6).
4. Once in place, tighten each end around the body for a comfortable and snug fit.
5. To remove pad, detach the hook fasteners tabs by pulling up on one end of the pad. Unwrap pad from body.

**NOTE:** Some wear of the exterior pad fabric may be seen over time; this is normal and will not affect the performance of the heating pad. Use caution when fastening the pad to prevent hook fasteners from touching clothing.

## HOW TO USE YOUR HEATING PAD

1. After controller is connected to pad, plug into 120V AC electrical outlet. Note: Controller must be connected to heating pad before plugging unit into outlet.
2. To turn the pad on, press the Control Button (Figure 4, letter A).
3. Select from 3 heat settings: low, medium, or high, by repeatedly pressing the control button until the desired heat setting indicator lights up (Figure 4, letter B). It is recommended to start on the low heat setting. Gradually move to the medium and high heat settings in 15 minute increments if more heat is desired. If at any point during use, you feel discomfort, select a lower heat setting or turn the unit to the off setting. Check skin under pad frequently while using.
4. To turn the pad off, press the Control Button (Figure 4, letter A) repeatedly until all of the heat setting indicator lights (Figure 4, letter B) turn off.

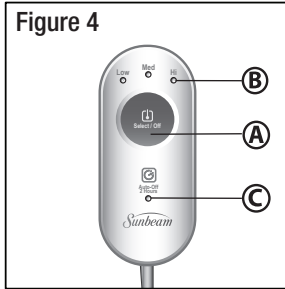- A. Control Button  
B. Heat Setting Indicator  
C. Auto-Off Indicator

### CAUTION:

- DO NOT SIT ON, LEAN AGAINST, CRUSH OR FOLD PAD DURING USE.
- DO NOT TO OVERLAP THE PAD. THE RED EDGES OF THE PAD SHOULD NEVER TOUCH THE PAD (Figure 1, number 6).
- DO NOT USE WHILE SLEEPING.
- SEE ADDITIONAL WARNINGS IN THE IMPORTANT SAFETY INSTRUCTIONS SECTION.

## AUTO-OFF FEATURE

Your heating pad controller includes an Auto-Off feature. This feature is designed to conserve energy by automatically turning off after 2 hours of continuous use. The Auto-Off Indicator (Figure 4, letter C) will flash when the auto-off feature is activated and the pad shuts off. To reset after Auto-Off, press the control button until the desired heat setting is reached.

**NOTE:** After a power outage, the Auto-Off feature automatically resets the time.

## HOW TO USE YOUR HEATING PAD WITH MOIST HEAT THERAPY

In some cases, you may find that moist heat therapy provides more comfort and relief. When using moist heat therapy, you should protect furniture and clothing from becoming damp.

1. Dampen the interior pad fabric (Figure 1, number 1) of the heating pad with water by patting the surface with a wet sponge or cloth, or use a spray bottle to apply a fine mist. Do not thoroughly soak pad.
2. Do not spray or dampen the pad connector (Figure 5).

**DANGER:** DO NOT IMMERSE THE PAD OR WIRING IN WATER OR PLACE UNDER RUNNING WATER. FAILURE TO COMPLY MAY RESULT IN PRODUCT FAILURE AND/OR PERSONAL INJURY. THE ELECTRICAL PAD CONNECTOR MUST BE THOROUGHLY DRY BEFORE ESTABLISHING ELECTRICAL CONNECTION.

Figure 5

## WASHING INSTRUCTIONS

DO NOT DRY CLEAN! DO NOT USE BLEACH! DO NOT USE WRINGER!

1. Disconnect pad from outlet. DO NOT PULL BY THE CORD.
2. Disconnect controller cord from pad by holding the pad side connector with one hand and pulling firmly on the cord side connector with the other hand.
3. Fold pad ends and attach hook fasteners to the exterior pad fabric before placing in the washing machine, or separate in a laundry wash bag.
4. Machine wash pad in cold water on gentle cycle or hand wash. Tumble dry on low heat and remove promptly from dryer or drape to dry.
5. DO NOT plug in pad and pad connector until completely dry.
6. For information on service, visit our website at [www.sunbeam.com](http://www.sunbeam.com).

## TROUBLE SHOOTING

Problem	Solution
No Heat / No Power	<ul style="list-style-type: none"> <li>• Verify tight connection between pad and controller. If disconnected, to restart unit, press the Control Button until the desired heat setting is reached.</li> <li>• Verify AC outlet is functional.</li> <li>• Verify no physical damage to the cord, controller or pad. If there is damage, discontinue use.</li> <li>• If failure continues, discontinue use. See warranty information to contact Sunbeam Consumer Service.</li> </ul>
Pad has No Heat & Controller 'Auto-Off' Indicator is Flashing	<ul style="list-style-type: none"> <li>• Auto-Off feature has been activated and set to turn the unit off after 2 hours of continuous use.</li> <li>• To re-start unit, press the Control Button until the desired heat setting is reached.</li> </ul>
Pad has a Hot Spot	<ul style="list-style-type: none"> <li>• Ensure the pad is not folded, bunched, overlapped or crushed in any area.</li> </ul>
'High' Heat Setting Indicator Light is Flashing	<p>The 'High' Heat Setting Indicator light will start flashing when a fault is detected. Follow these steps to reset:</p> <ol style="list-style-type: none"> <li>1. Unplug the controller from the AC outlet.</li> <li>2. Verify tight connection between pad and controller.</li> <li>3. After 10 seconds, plug the controller in the AC outlet.</li> <li>4. Press the Control Button until the desired heat setting is reached.</li> </ol> <p>If flashing light persists, an error has occurred. Discontinue use of the product and see warranty information to contact Sunbeam Customer Service.</p>

## FCC CERTIFIED

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

**Note:** This product has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This product generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the product and receiver.
- Connect the product into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

## 2-YEAR LIMITED WARRANTY

Sunbeam Products, Inc. ("Sunbeam") or if in Canada, Newell Industries Canada ULC ("NIC") warrants that for a period of two (2) years from the date of initial retail purchase, this product will be free from defects in material and workmanship, when used in accordance with the instructions provided.

Sunbeam (or NIC), at its option, will repair or replace this product or any component of the product found to be defective during the warranty period.

Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. Do NOT attempt to repair or adjust any electrical or mechanical functions on this product. Doing so will void this warranty.

This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. Sunbeam (or NIC) dealers, service centers, or retail stores selling Sunbeam (or NIC) products do not have the right to alter, modify or any way change the terms and conditions of this warranty.

This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than Sunbeam (or NIC) or an authorized Sunbeam (or NIC) service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

### What are the limits on Sunbeam's (or NIC's) Liability?

Sunbeam (or NIC) shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition. Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty. Sunbeam (or NIC) disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise. Sunbeam (or NIC) shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party. Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you. This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

### How to Obtain Warranty Service

#### In the U.S.A.

If you have any question regarding this warranty or would like to obtain warranty service, please call 1 800 435-1250 and a convenient service center address will be provided to you.

#### In Canada

If you have any question regarding this warranty or would like to obtain warranty service, please call 1 800 667- 8623 and a convenient service center address will be provided to you.

In the U.S.A., this warranty is offered by Sunbeam Products, Inc. located in Boca Raton, Florida 33431.

In Canada, this warranty is offered by Newell Industries Canada ULC, located at 20B Hereford Street, Brampton, Ontario L6Y 0M1. If you have any other problem or claim in connection with this product, please write our Consumer Service Department.

**PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE.**

For inquiries regarding recycling and proper disposal of this product, please contact your local waste management facility.

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