
Health o meter®

WEIGHT TREND TRACKING SCALE

User's Manual
Model HDM037

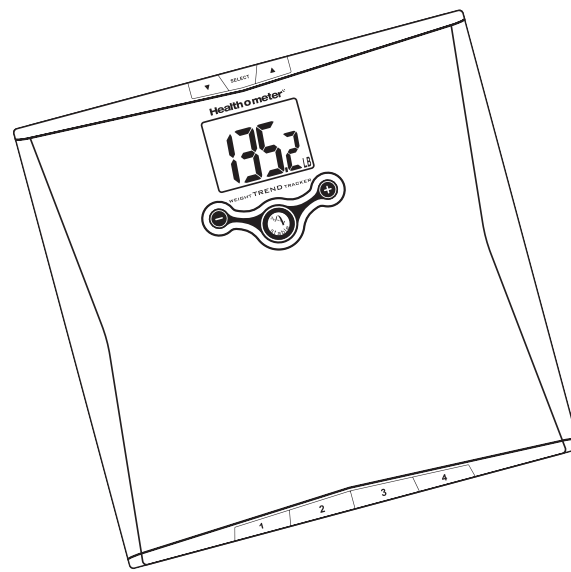www.healthometer.com

WELCOME!

Congratulations on your purchase of the Health o meter® wellness monitoring product. Not only will your scale measure your weight with accuracy every time, but it also includes an exclusive Trend Tracker feature that is a **quick and easy** visual indication of how your weight is trending. It also includes up to four user profiles.

ICONS AND INDICATORS

- Body Mass Index (BMI) Indicator
- User and Reset Icons
- GOAL Selection: LOSE or GAIN
- Step On/Off Indicators
- Unit of Height Measurement: Imperial System - Feet (") and Inches (") or Metric System - Centimeters (CM)
- Unit of Weight Measurement: Pounds (LB) or Kilograms (KG)
- Weight and BMI Readout
- Weight Trend Tracker Indicator-GAINED
- Weight Trend Tracker Indicator-STABLE
- Weight Trend Tracker Indicator-LOST


SELECTING LB/KG AND RESETTING SCALE

The unit is automatically set up to measure in pounds. If you want to switch the unit of measurement then follow the steps below.

- Press and hold the **SELECT** button for a few seconds when the scale is off.

- Use the **↑** and **↓** to scroll through **LB** or **KG**. Press **SELECT** to confirm.

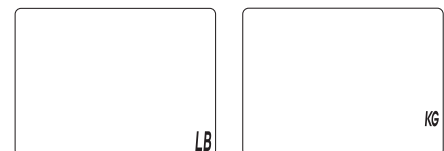
- Use the **↑** and **↓** to scroll through **YES** or **NO** to **RESET**. Press **SELECT** to confirm.


CAUTION Confirming to reset will erase **ALL** user information.

FIRST TIME USER SET UP

- Press the desired **USER** button and the **HEIGHT** will begin blinking.
- Press the **↑** or **↓** button to choose the desired height in **FEET** or **Centimeters (CM)**. Then press the **SELECT** button.
- If using the Imperial unit of measurement (feet/inches), press the **↑** or **↓** button to choose the desired **INCHES**. Then press the **SELECT** button.
- If you want to select a **GOAL**, press the **↑** or **↓** button to select **YES**. Then press the **SELECT** button.

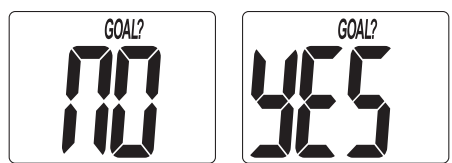- If you selected **YES** in step 4, press the **↑** or **↓** button to select **GAIN** or **LOSE**. Then press the **SELECT** button.

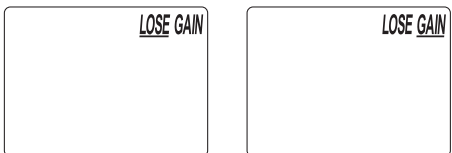- If you selected **NO** or did not finish step 4, the Trend Tracker feature will compare the weight trend to the user's average weight.

The user set up is complete. The scale will calibrate for a few seconds and then be ready for use.

NOTE: To reset a **specific user**, please follow the next steps.

- Press the desired user number button.
- When user number is displayed, immediately press and hold **SELECT** button.
- Press **↑** and **↓** button to choose **YES** or **NO**. Then press **SELECT** to confirm.

TREND TRACKER FEATURE

The trend tracker feature provides the user with a quick indication of how their weight is trending.

Green LED: will indicate if a user's weight is **in line** with the desired trend (GAIN or LOSE)

Red LED: will indicate if a user's weight is **not in line** with the desired trend (GAIN or LOSE)

Amber LED: will indicate if a user's weight is **stable**.


	Lost	Stable	Gained
For users weighing less than 100 LB (45 Kg)	Your current weight is 2 lbs. ≤ below your average weight.	Your current weight is within 2 lbs. of your average weight.	Your current weight is 2 lbs. ≥ above your average weight.
For users weighing 100 LB (45 Kg) or more	Your current weight is 2% ≤ below your average weight.	Your current weight is within 2% of your average weight.	Your current weight is 2% ≥ above your average weight.

MONITORING YOUR WEIGHT

- Press one of the **USER** buttons. The display will show the user number followed by height. The scale will calibrate itself and will show **0.0** and the unit of measurement (LB or KG).

- Step on the scale and position your feet in a manner that will distribute your weight over the entire scale platform. Stand still while reading the display. The display will show your weight, BMI and the trend light.

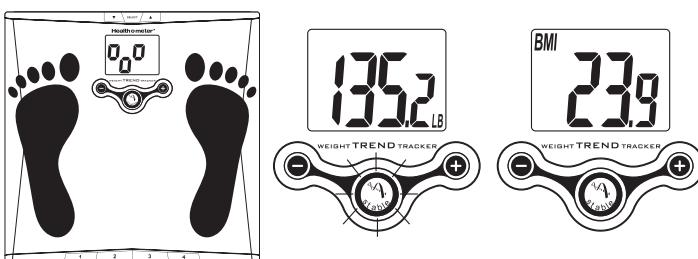
NOTE: In order to preserve battery life, your scale will automatically turn off after a short period of inactivity.

WEIGHING WITHOUT SELECTING A USER

CAUTION By not selecting a user, the Trend Tracker feature will not be activated.

- Place the scale on a flat surface. A hard surface is best for maximum accuracy.
- Step onto the scale and position your feet so that your weight is distributed evenly over the scale platform. The LCD display will flash "0.0" several times for a few seconds before displaying your weight. Stand still while reading the display.

NOTE: This scale offers both "Instant On" (Step onto scale) or "Tap On" (Tap on lower right corner and wait for "0.0" and "Step On"). During the "Instant On" weighing, if you see "E", the scale has recalibrated. Please step on the scale to weigh yourself again.

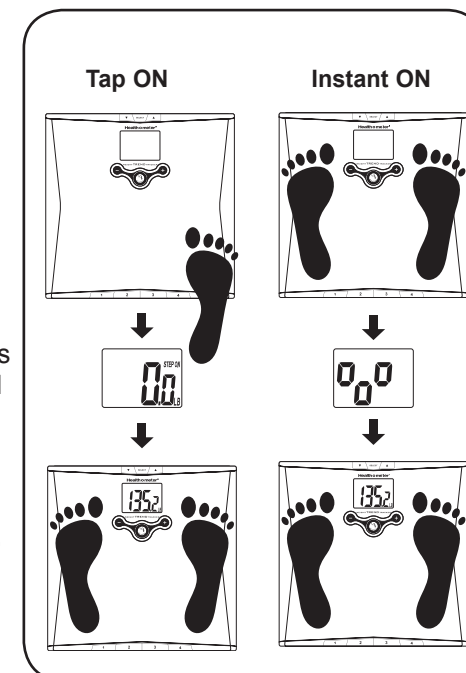
REPLACING THE BATTERIES

This unit includes four (4) AA batteries. When "Lo" appears on the display, replace all batteries.

CAUTION Do not mix old and new batteries. Do not mix alkaline, standard (Carbon-Zinc) or rechargeable (Nickel-Cadmium) batteries. Removing the battery will not affect the recorded data and user profile.

Batteries Replacement

- Locate the battery cover on the back of the scale. Push the locking tab and pull the battery cover up.
- Use the tip of a pen to loosen each battery. Hold the scale at an angle to allow the battery to slide out.
- Insert four new batteries, pressing them firmly into their slot.
- Replace the battery cover.

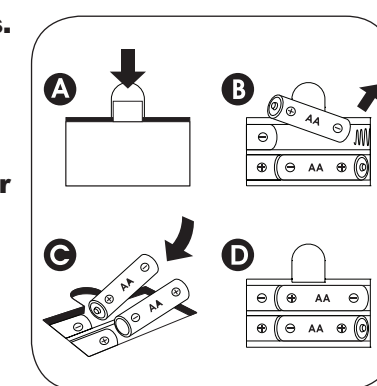CARING FOR YOUR SCALE


Periodically wipe the scale surface down with a DRY cloth to remove dust.

CAUTION Weight exceeding the capacity (400 LB/181 kg) may damage your scale. This product is designed for personal use only. Not for commercial or industrial use. This unit is not waterproof; avoid contact with excessive moisture.

TROUBLESHOOTING

If "E" appears on the display:

- The scale cannot detect a stable weight. Stand still while weighing.
- The weight has exceeded the scale's capacity (400 lb / 181 kg).

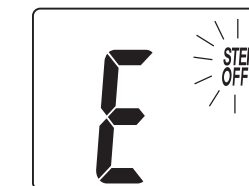FCC STATEMENT

POTENTIAL FOR RADIO/TELEVISION INTERFERENCE

This product has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

The product generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that the interference will not occur in a particular installation. If the product does cause harmful interference to radio or television reception, which can be determined by turning the product on or off, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna;
- Increase the separation between the product and the receiver;
- Connect the product into an outlet on a circuit different from that to which the receiver is connected;
- Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

5 YEAR LIMITED WARRANTY

PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE

Sunbeam Products, Inc. doing business as Jarden Consumer Solutions or if in Canada, Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions (collectively "JCS") warrants that for a period of five years from the date of purchase, this product will be free from defects in material and workmanship. JCS, at its option, will repair or replace this product or any component of the product found to be defective during the warranty period. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. Do NOT attempt to repair or adjust any electrical or mechanical functions on this product. Doing so will void this warranty. This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. JCS dealers, service centers, or retail stores selling JCS products do not have the right to alter, modify or in any way change the terms and conditions of this warranty. This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than JCS or an authorized JCS service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

What are the limits on JCS's Liability?

JCS shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition. Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty. JCS disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise. JCS shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party. Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you. This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

How to Obtain Warranty Service

In the U.S.A.: If you have any questions regarding this warranty or would like to obtain warranty service, please call 1 800-672-5625 and a convenient service center address will be provided to you.

In Canada: If you have any questions regarding this warranty or would like to obtain warranty service, please call 1 800-667-8623 and a convenient service center address will be provided to you.

In the U.S.A., this warranty is offered by Sunbeam Products, Inc. doing business as Jarden Consumer Solutions located in Boca Raton, Florida 33431. In Canada, this warranty is offered by Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions, located at 20 B Hereford Street, Brampton, Ontario L6Y 0M1. If you have any other problem or claim in connection with this product, please write our Consumer Service Department.

PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE

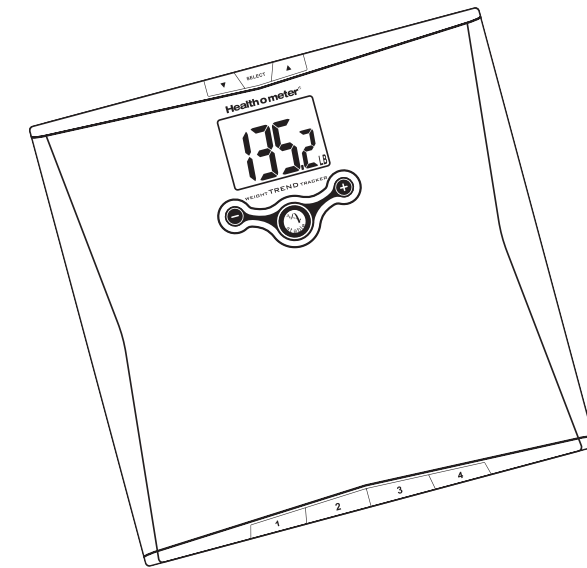
© 2010 Sunbeam Products, Inc. doing business as Jarden Consumer Solutions. All rights reserved. Distributed by Sunbeam Products, Inc. doing business as Jarden Consumer Solutions, Boca Raton, Florida 33431.

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