Renue™ Back & Body Pad

Model: 889


IMPORTANT SAFETY INSTRUCTIONS – DO NOT DESTROY

WARNINGS AND SAFEGUARDS

AS WITH ALL ELECTRICAL PRODUCTS, MISUSE OF PRODUCT OR FAILURE TO PROPERLY FOLLOW THE INSTRUCTIONS MAY CAUSE OVERHEATING, FIRE OR PERSONAL INJURY. PLEASE READ THE PRODUCT LABEL AND ALL OF THE INSTRUCTIONS BEFORE USING YOUR HEATED PRODUCT.

1. DO NOT USE THIS HEATED PRODUCT WITH AN INFANT, A CHILD, AN INCAPACITATED OR IMMOBILE PERSON, A PARAPLEGIC, A QUADRIPLAGIC, A DIABETIC, ONE INSENSITIVE TO HEAT SUCH AS A PERSON WITH POOR BLOOD CIRCULATION, OR ANYONE WHO CANNOT CLEARLY UNDERSTAND THE INSTRUCTIONS AND/OR OPERATE THE CONTROLS.
2. WHEN USED ON A BED, DO NOT ROUTE ELECTRICAL CORDS BETWEEN THE MATTRESS AND BOX SPRING. DO NOT ALLOW THE CORDS TO BE PINCHED.
3. DO NOT USE WITH DAMAGED CORDS. DAMAGED CORDS CAN LEAD TO FIRE. DO NOT TUCK IN, TRAP, FOLD, CROSS OR PINCH CORDS BETWEEN CHAIR'S SEAT AND BACK OR ARMS, AGAINST WALLS, BETWEEN BED SLATS, MATTRESS OR SPRINGS, FOOTBOARDS OR BED FRAME.
4. DO NOT DRY CLEAN THIS HEATED PRODUCT. CLEANING LIQUID AND BLEACH MAY DAMAGE THE INSULATION OF THE HEATING ELEMENT. IMPROPER LAUNDERING MAY CAUSE PROBLEMS WITH THE HEATING WIRE WHICH CAN LEAD TO FIRE.
5. DO NOT TUCK IN, PINCH, FLEX OR EXPOSE THE HEATING WIRE TO RUBBING. A DAMAGED HEATING WIRE MAY RESULT IN POSSIBLE OVERHEATING OR FIRE. DO NOT TUCK THE HEATED PRODUCT OR SUPPLY CORD IN A SOFA, ARMCHAIR OR OTHER FURNITURE.
6. TURN OFF OR UNPLUG THE CONTROL FOR THE HEATED PRODUCT WHEN NOT USING.
7. WRAPPING THE CORD AROUND THE CONTROLLER MAY DAMAGE THE CORD. LOOP LOOSELY WHEN STORING.
8. CONTROLLER IS NOT WATERPROOF. DO NOT IMMERSE IN WATER.

9. KEEP THE CONTROL AWAY FROM DRAFTY AREAS, SUCH AS BY AN OPEN WINDOW, WHICH MAY CAUSE IT TO BECOME DAMP OR WET.
10. DO NOT USE PINS; THEY MAY DAMAGE THE ELECTRIC WIRING.
11. BE SURE THAT THIS HEATED PRODUCT IS USED ONLY ON AN AC SUPPLY CURRENT CIRCUIT OF 110-120 VOLTS. DO NOT USE WITH AN EXTENSION CORD.
12. THE PRODUCT MAY EXPERIENCE INCREASED HEAT IN AREAS WHERE HEATED PRODUCT IS BUNCHED OR FOLDED WHILE IN USE.
13. KEEP ALL DOGS, CATS AND OTHER PETS AWAY FROM THIS HEATED PRODUCT. THIS IS AN ELECTRICAL APPLIANCE THAT THEY CAN DAMAGE OR BE INJURED BY.
14. DO NOT USE THIS HEATED PRODUCT ON A WATERBED, BUNK BED, MECHANICALLY ADJUSTABLE BED OR A CHILD'S CRIB/BED.
15. DO NOT USE WITH ROCKING CHAIRS, RECLINERS, CHAIRS WITH CASTERS OR WHEELS, OR OTHER FURNITURE SUBJECT TO MOVEMENT. DO NOT WRAP PRODUCT AROUND CONTROL.
16. DO NOT USE ANOTHER ELECTRONIC HEATING DEVICE IN CONJUNCTION WITH YOUR HEATED PRODUCT.
17. IF THESE IMPORTANT INSTRUCTIONS ARE NOT FOLLOWED, THE HEATED PRODUCT WILL NOT OPERATE PROPERLY.
18. NOT FOR HOTEL, MOTEL OR OTHER INSTITUTIONAL USE.
19. DO NOT USE WITH GENERATOR POWER.
20. IF YOU SEE THE HEATED PRODUCT OPERATING IMPROPERLY, UNPLUG IT IMMEDIATELY AND CALL SUNBEAM CONSUMER CARE IN THE U.S AT 800-435-1250 OR IN CANADA AT 800-667-8623.

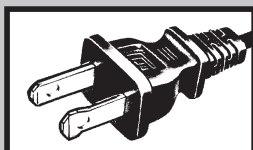THIS APPLIANCE HAS A POLARIZED PLUG (ONE BLADE IS WIDER THAN THE OTHER). TO REDUCE THE RISK OF ELECTRIC SHOCK, THIS PLUG IS INTENDED TO FIT IN A POLARIZED OUTLET ONLY ONE WAY. IF THE PLUG DOES NOT FIT FULLY IN THE OUTLET, REVERSE THE PLUG. IF IT STILL DOES NOT FIT, CONTACT A QUALIFIED ELECTRICIAN TO INSTALL THE PROPER OUTLET.

DO NOT ATTEMPT TO MODIFY THIS PLUG OR DEFEAT THIS SAFETY FEATURE IN ANY WAY.

THERAPEUTIC LOW HEAT

Unlike traditional heating pads, your Renue™ Back & Body Warming Pad is designed to provide a low-level therapeutic heat over an extended period of time. When used as intended – on a chair or on a bed – it will give off consistent and soothing low heat, providing therapeutic relief from aches and pains. The heat this pad generates can best be experienced by seating or lying on top of the product.

SET-UP & USAGE INSTRUCTIONS

Your Renue™ Back & Body Warming Pad is designed to be used on a chair or on bed for a continuous heat therapy. It provides relief from aches and pains by giving off consistent low level, soothing heat while you work or relax. This pad can also be used as a proactive heat therapy – to reduce stress and tension before it hurts!

FOR USE ON A CHAIR:

1. Place the Renue™ Back & Body Pad on a chair, with softer fleece fabric facing out. The cord connector should be on the right side of the chair's seat, and the two loop holes for straps should be at the top corners of the chair's back (Figure 1).
2. Secure the pad on the chair with provided flexible Velcro straps (Figures 2a-b). Depending on personal preference and the size and type of chair, the orientation of the straps may change. Place the straps in a way that allows the pad to be attached securely to the chair while still being comfortable to sit on.
3. As an optional feature, your Renue™ Back & Body pad can be rolled up for lumbar support, allowing for added comfort throughout the day. Just roll up the seat portion of the pad and secure with the Velcro strip on the underside of the pad (Figure 3).
4. Make sure Control Connector is securely attached to the Pad Connector (Figure 4). Push both parts together until they click for a tight fit.
NOTE: Control must be connected to heating pad before plugging unit into outlet.
5. After control is connected to pad, plug into 120V AC electrical outlet. For convenience, the pad features a snap that can be used to keep cord & controller handy, but out of the way (Figure 5).
NOTE: To avoid any tangling or pinching, do not tuck in the cord between chair's seat and back or arms.
6. To turn the pad on, press the Power Button. Then select from 4 heat settings by scrolling to the desired setting with the Power Button from "warm" to "high" (Figure 6).
7. To turn the pad OFF, scroll through the remaining heat settings by pressing the Power Button sequentially until no lights remain lit.

FOR USE ON A BED:

1. The Renue™ Back & Body pad can be placed on top of or underneath a fitted sheet.
2. Place the pad on a bed in the position that best allows heat to target the affected area of aches or pains. The softer fleece fabric should be facing up (Figure 7).
3. Make sure Control Connector is securely attached to the Pad Connector (Figure 4). Push both parts together until they click for a tight fit.
NOTE: Control must be connected to heating pad before plugging unit into outlet.
4. After control is connected to pad, plug into 120V AC electrical outlet.
5. Run the control cord on the floor under the bed. To avoid any tangling or pinching, do not sandwich cord between the mattress and box spring. For your convenience, the pad features a snap that can be used to keep the cord & controller handy, but out of the way.
6. To turn the pad on, press the Power Button. Then select from 4 heat settings by scrolling to the desired setting with the Power Button from "warm" to "high" (Figure 6).
7. To turn the pad OFF, scroll through the remaining heat settings by pressing the Power Button sequentially until no lights remain lit.

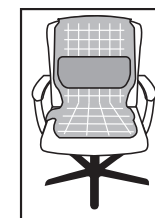
Figure 1

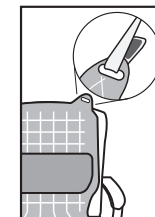
Figure 2a

Figure 2b

Figure 3

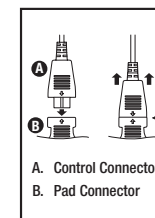
Figure 4

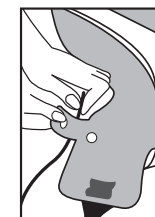
Figure 5

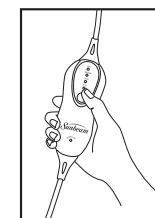
Figure 6

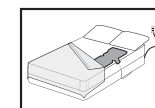
Figure 7

AUTO-OFF FEATURE

Your pad controller includes an Auto-Off feature. This feature is designed to conserve energy by automatically turning off after 8 hours of continuous use.

NOTE: After a power outage, the Auto-Off feature automatically resets the timer.

CARE AND CLEANING

BEFORE WASHING OR STORING YOUR HEATED PRODUCT:

1. Disconnect power cord from wall BEFORE detaching control cord from the heated product.
 2. Ensure that no cords are still connected to the pad and none of the wires have worn through the fabric of the pad.
 3. Carefully follow all washing instructions or your Renue™ Back & Body pad may shrink and/or cause a fire hazard.
- DO NOT submerge control in water.
 - DO NOT dry clean your heated product. Dry cleaning solvents may damage the wires.
 - DO NOT use bleach or cleaning fluids.
 - DO NOT use a wringer.
 - DO NOT iron your heated product.
 - DO NOT use mothballs or sprays of any kind.
 - DO NOT reconnect the heated product or use it until it is COMPLETELY dry.
 - DO NOT wash or put into water the power cord or any of the product's cords.

Taking care of your Renue™ Back & Body Pad is easy as long as you follow a few simple steps.

MACHINE WASHING

1. Presoak the pad for 15 minutes in mild soap and lukewarm water.
2. Wash in mild soap and lukewarm water on the "delicate" or "gentle" cycle in the washer for two minutes. Remember, do not use bleach.
3. Rinse in fresh, lukewarm water.
4. Spin dry in the washing machine. Then gently stretch to its original size.

MACHINE DRYING

CAUTION: Do not use commercial dryers or the dryers at your local laundromat. They get too hot and will damage the wiring of the heated product.

1. Preheat the dryer on the "low" temperature setting for 2 minutes.
2. Place pad into the dryer. Remove the pad while it is still damp, then gently stretch pad to its original size.
3. Drape pad over a clothesline or a shower rod. DO NOT use clothespins.
4. Pad must be COMPLETELY dry before reconnecting the power plug or control plug.

HAND WASHING

1. Presoak the pad for 15 minutes in mild soap and lukewarm water.
2. Gently squeeze the soap suds through a few times.
3. Rinse in fresh, lukewarm water.

HAND DRYING

1. Gently stretch the pad back to its original size.
2. Drape over a clothesline or hang over a shower rod. DO NOT use clothespins.
3. Pad must be COMPLETELY dry before reconnecting the power plug or control plug.

5-YEAR LIMITED WARRANTY

Sunbeam Products, Inc. doing business as Jarden Consumer Solutions or if in Canada, Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions (collectively "JCS") warrants that for a period of five years from the date of purchase, this product will be free from defects in material and workmanship. JCS, at its option, will repair or replace this product or any component of the product found to be defective during the warranty period. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. Do NOT attempt to repair or adjust any electrical or mechanical functions on this product. Doing so will void this warranty. This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. JCS dealers, service centers, or retail stores selling JCS products do not have the right to alter, modify or any way change the terms and conditions of this warranty. This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than JCS or an authorized JCS service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

What are the limits on JCS's Liability?
JCS shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition. Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty.

JCS disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise. JCS shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party. Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you. This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

How to Obtain Warranty Service
In the U.S.A. - If you have any question regarding this warranty or would like to obtain warranty service, please call 1-800-435-1250 and a convenient service center address will be provided to you.

In Canada - If you have any question regarding this warranty or would like to obtain warranty service, please call 1-800-667-8623 and a convenient service center address will be provided to you.

In the U.S.A., this warranty is offered by Sunbeam Products, Inc. doing business as Jarden Consumer Solutions located in Boca Raton, Florida 33431. **In Canada**, this warranty is offered by Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions, located at 20 B Herford Street, Brampton, Ontario L6Y 0M1. If you have any other problem or claim in connection with this product, please write our Consumer Service Department.

PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE.

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©2012 Sunbeam Products, Inc. operando bajo el nombre de Jarden Consumer Solutions. Todos los derechos reservados. Distribuido por Sunbeam Products, Inc. operando bajo el nombre de Jarden Consumer Solutions, Boca Raton, Florida 33431. Impreso en China. DW-080112 PN: 156796