Health o meter™

DIGITAL TAPE MEASURE

User's Manual
Model HDTM012-69

www.healthometer.com

WELCOME!

Congratulations on your purchase of this Health o meter™ digital tape measure. This product will help to measure and store up to eight body part circumferences.

PLEASE READ CAREFULLY AND SAVE

For information on service or warranty, contact Sunbeam Products, Inc. at 1-800-672-5625 or go to www.healthometer.com.

WARNING

Not a toy.

Not to be used by children without adult supervision.


This product is designed for personal use to measure body circumferences only. Not for commercial or industrial use, or linear measurement of length.

Do not retract the tape too fast.

ABOUT YOUR HEALTH O METER TAPE MEASURE

The Health o meter™ Digital Tape Measure helps you to measure and record up to eight body circumferences.

Keys and Icons

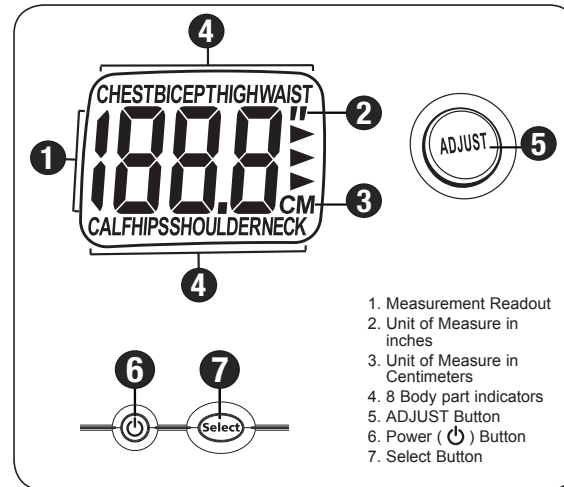
OPERATING INSTRUCTIONS

1. Press the **POWER** (⏻) button to turn the unit on. Make sure the tape is fully retracted.
2. The screen will automatically prompt a body part to be measured.
3. Pull the tape and wrap it around the body part.
4. Insert the locking pin into the tape receptacle.
5. Push and hold the **ADJUST** button until the tape has fully retracted around the body part.
6. Press "**Select**" to confirm the measurement. The unit will make an audible beep.
7. Press "**Select**" to save the reading to memory and move to the next body part. The unit will make an audible beep.
8. Fully retract tape. This will prepare the unit for the next measurement and ensure accuracy.
9. Repeat the same set of steps for each of the body parts.

NOTE: if you don't want to measure a certain body part then press the "Select" button to move to the next body part.

Measuring sites:

- Neck:** While standing, measure neck at largest girth above the Adam's apple
- Shoulder:** While standing, measure the greatest circumference around your shoulders. Take measurement around the widest part of your arm. Arm can be bent or straight.
- Bicep:** While standing, measure the circumference of your chest just above the nipple area. Make sure you have exhaled prior to taking measurement.
- Chest:** While standing, measure at the midway point between the top of your hip and the bottom of your rib cage.
- Hips:** While standing, measure at the largest girth around your hips and buttocks.
- Thigh:** While standing, measure the widest girth below the buttocks.
- Calf:** While seated, measure at its widest girth.

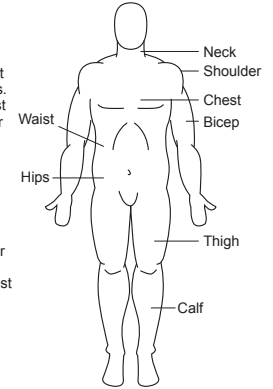Directions:

- It is recommended to take measurements in a relaxed state. Do not flex while taking the measurements.
- Measurements should be taken prior to any physical activity.
- Measurements should be taken on the right side of the body.

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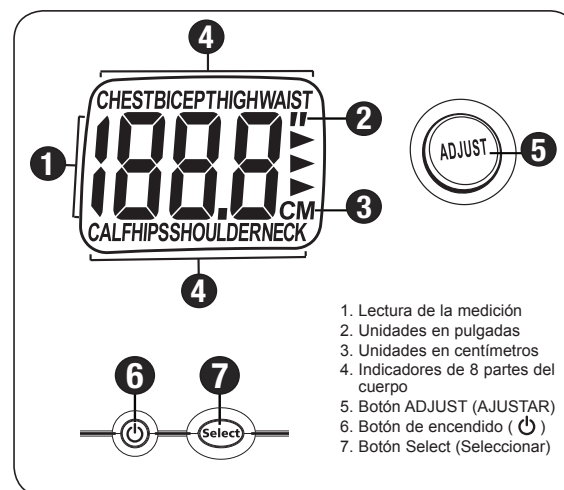
CAUTION This product is designed for personal use only. Not for commercial or industrial use. This unit is not waterproof; avoid contact with excessive moisture.

Recording Measurements

In order to save the measurement for a body part, press the "Select" button while the displayed reading is blinking. The measurement will be displayed on the screen when you are ready to measure that body part again.

NOTE: All saved readings will be reset to zero when battery is replaced.

NOTE: The previously stored measurement will be replaced when a new measurement is saved.

Viewing Prior Measurements

To scroll through prior measurements press the "Select" button to cycle through the different body parts.

Calibrating the Unit

If you want to zero out the unit press and hold the POWER () and "Select" buttons. Make sure that the tape is fully retracted when you are doing this.

Changing Units of Measurement

You can switch between inches and centimeters by pressing and holding the "Select" button.

Turning the Unit Off

To turn the unit off, press and hold the POWER () button.

TROUBLESHOOTING

Error Display ("E")

If the readout is a negative value, "E" will appear on the display:

- To reset the unit, turn the unit off by pressing and holding the POWER () button. Make sure that the tape is fully retracted. Follow the OPERATING INSTRUCTIONS to re-start the measurements again.

- To zero out the unit, follow the above section of "Calibrating the Unit".

NOTE: Display will turn off automatically in 2 minutes if no button is pressed.

Battery Low Display ("Lo")


If the battery power is weak, "Lo" will appear on the display, follow the subsequent section of "REPLACING THE BATTERY" to replace with new battery.

REPLACING THE BATTERY

- Locate the battery cover on the back of the tape measure.
- Use a flat head screwdriver or coin to rotate the battery cover counter clockwise to the Unlock () position.
- Remove the used battery and insert a new Lithium battery (CR2032-3V).
- Rotate the battery cover clockwise to the Lock () position.

NOTE: The Lithium coin battery(s) used in this product contains perchlorate material. Special handling may apply in California. Go to: www.dtsc.ca.gov/hazardouswaste/perchlorate for more information.

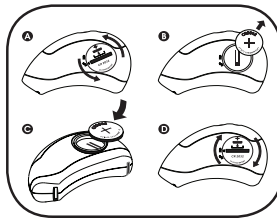CARING FOR YOUR TAPE MEASURE

Periodically wipe the surface of the tape measure with a DRY cloth to remove dust.

POTENTIAL FOR RADIO/TELEVISION INTERFERENCE


This product has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

The product generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that the interference will not occur in a particular installation. If the product does cause harmful interference to radio or television reception, which can be determined by turning the product on or off, the user is encouraged to try to correct the interference by one or more of the following measures: a) Reorient or relocate the receiving antenna; b) Increase the separation between the product and the receiver; c) Connect the product into an outlet on a circuit different from that to which the receiver is connected; d) Consult the dealer or an experienced radio/TV technician for help. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

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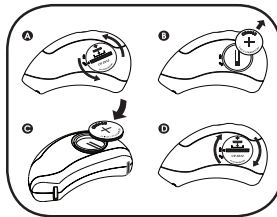5 YEAR LIMITED WARRANTY

PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE.

Sunbeam Products, Inc. doing business as Jarden Consumer Solutions or if in Canada, Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions (collectively "JCS") warrants that for a period of five years from the date of purchase, this product will be free from defects in material and workmanship. JCS, at its option, will repair or replace this product or any component of the product found to be defective during the warranty period. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. DO NOT attempt to repair or adjust any electrical or mechanical functions on this product. Doing so will void this warranty. This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. JCS dealers, service centers, or retail stores selling JCS products do not have the right to alter, modify or in any way change the terms and conditions of this warranty. This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than JCS or an authorized JCS service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

What are the limits on JCS's Liability?

JCS shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition. Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty. JCS disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise. JCS shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party. Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you. This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

How to Obtain Warranty Service

In the U.S.A.: If you have any questions regarding this warranty or would like to obtain warranty service, please call 1 800-672-5625 and a convenient service center address will be provided to you.

In Canada: If you have any questions regarding this warranty or would like to obtain warranty service, please call 1 800-667-8623 and a convenient service center address will be provided to you.

In the U.S.A., this warranty is offered by Sunbeam Products, Inc. doing business as Jarden Consumer Solutions located in Boca Raton, Florida 33431. In Canada, this warranty is offered by Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions, located at 20 B Hereford Street, Brampton, Ontario L6Y 0M1. If you have any other problem or claim in connection with this product, please write our Consumer Service Department.

PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE.

© 2012 Sunbeam Products, Inc. doing business as Jarden Consumer Solutions. All rights reserved. Distributed by Sunbeam Products, Inc. doing business as Jarden Consumer Solutions, Boca Raton, Florida 33431. Printed in China.

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