Sunbeam

GoHeat™ Cordless Heating Pad

Models: 2000-825 & 2003-825

IMPORTANT SAFETY INSTRUCTIONS DO NOT DESTROY

DANGER: TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK AND FIRE, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:


- READ ALL INSTRUCTIONS.
- DO NOT USE WHILE SLEEPING OR ON SOMEONE WHO IS SLEEPING.
- DO NOT USE ON AN INFANT OR ON AN ANIMAL.
- THIS PAD IS NOT TO BE USED ON OR BY AN INVALID, SLEEPING, HELPLESS OR UNCONSCIOUS PERSON, A PERSON WITH DIABETES, A PERSON WITH POOR BLOOD CIRCULATION OR INSENSITIVITY TO HEAT, A PERSON INCAPACITATED BY MEDICATION.
- THIS APPLIANCE IS NOT TO BE USED ON OR BY CHILDREN. CLOSE SUPERVISION IS NECESSARY WHEN ANY APPLIANCE IS USED NEAR CHILDREN. CHILDREN SHOULD NOT PLAY WITH THE APPLIANCE.
- DO NOT USE ON AREAS OF INSENSITIVE SKIN.
- BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING. CHECK SKIN UNDER PAD FREQUENTLY.
- DO NOT USE IN AN OXYGEN ATMOSPHERE.
- DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.
- DO NOT WRAP PAD TOO TIGHT AROUND THE BODY. PAD SHOULD BE SNUG, BUT SHOULD NOT CUT OFF CIRCULATION.
- DO NOT WRAP PAD AROUND OR NEAR THE NECK TO AVOID PERSONAL INJURY.
- DO NOT SIT ON, AGAINST OR CRUSH PAD – AVOID SHARP FOLDS. PLACE PAD ON TOP OF AND NOT UNDER THE PART OF BODY NEEDING HEAT.
- UNPLUG WHEN NOT IN USE. NEVER LEAVE APPLIANCE UNATTENDED, ESPECIALLY IF CHILDREN ARE PRESENT.
- NEVER PULL THIS PAD BY THE SUPPLY CORD. DO NOT USE THE CORD AS A HANDLE. DO NOT PINCH CORD.
- DO NOT TUCK IN, TRAP, FOLD, CROSS, PINCH OR WRAP CONTROL OR CORD WITHIN HEATING PAD DURING USE.
- USE CAUTION WITH STRAP END TO PREVENT FASTENING TAB FROM DAMAGING CLOTHING OR OTHER PROPERTY.
- DO NOT PLACE ANYTHING OVER THE HEATING PAD WHILE IN USE OR PLUGGED IN, SUCH AS A TOWEL, BLANKET, CLOTHING, OR ANY OTHER INSULATING MATERIAL, AND NEVER ALLOW HEATING PAD TO WRAP AROUND ITSELF.
- CAREFULLY EXAMINE PAD, COVER, CORD, AND BATTERY UNIT/HEATING PAD CONTROLLER BEFORE EACH USE. DISCARD THE PAD IF ANY SIGN OF DETERIORATION, WEAR OR DAMAGE. DO NOT TAMPER WITH THIS PAD IN ANY WAY. THERE ARE NO USER SERVICABLE PARTS. IF PAD DOES NOT FUNCTION SATISFACTORILY, SEE WARRANTY FOR CUSTOMER CONTACT INFORMATION.
- LOOP CORD LOOSELY WHEN STORING, TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS. ALLOW TO COOL BEFORE STORING AND DO NOT CREASE PAD BY PLACING ITEMS ON TOP OF IT. STORE IN THE ORIGINAL PACKAGING, IN COOL, DARK, DRY PLACE.
- USE ONLY THE PROVIDED BATTERY UNIT/HEATING PAD CONTROLLER AND PROVIDED POWER SUPPLY (SUNBEAM® CLASS 3S1P CHARGER) WITH THIS HEATING PAD.
- WHEN NOT IN USE, DISCONNECT HEATING PAD FROM THE BATTERY UNIT/HEATING PAD CONTROLLER BEFORE MAKING ANY ADJUSTMENTS, CHANGING ACCESSORIES, OR STORING. SUCH PREVENTATIVE SAFETY MEASURES REDUCE THE RISK OF STARTING THE APPLIANCE ACCIDENTALLY.

- DO NOT EXPOSE BATTERY UNIT/HEATING PAD CONTROLLER OR APPLIANCE TO FIRE OR EXCESSIVE TEMPERATURE. EXPOSURE TO FIRE OR TEMPERATURE ABOVE 266°F (130°C) MAY CAUSE EXPLOSION.
- DO NOT USE THIS PAD WITH LINIMENT, SALVE OR OINTMENT PREPARATIONS. SKIN BURNS COULD RESULT.
- THIS APPLIANCE IS NOT INTENDED FOR MEDICAL USE IN HOSPITALS.
- SAVE THESE INSTRUCTIONS.

HOUSEHOLD USE ONLY

HEATING PAD COMPONENTS

Figure 1
Models: 2000-825 & 2003-825

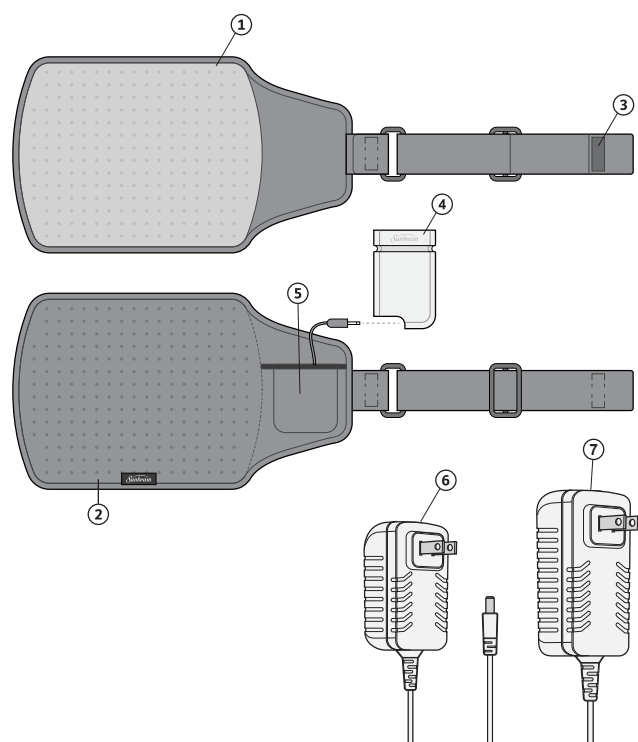
- Body Side (Lighter Colored Fabric)
- Exterior Side (Darker Colored Fabric)
- Fastening Tab
- Battery Unit/Heating Pad Controller (Also sold separately as Model: 2000-001)
- Integrated Pocket for Battery Unit/Heating Pad Controller
- Standard Charger
- Quick Charger (included with Model: 2003-825, also sold separately as Model: 2003-001)

WHERE TO USE YOUR HEATING PAD

The benefits of heat have been clinically proven to soothe sore muscles, reduce stiffness, and promote healing. The GoHeat™ Cordless Heating Pad provides the freedom to treat pain anytime, anywhere. The adjustable strap can be used to easily secure the pad around any area of pain. It is ideal for use on the back, arms, legs, and shoulders.

Figure 2

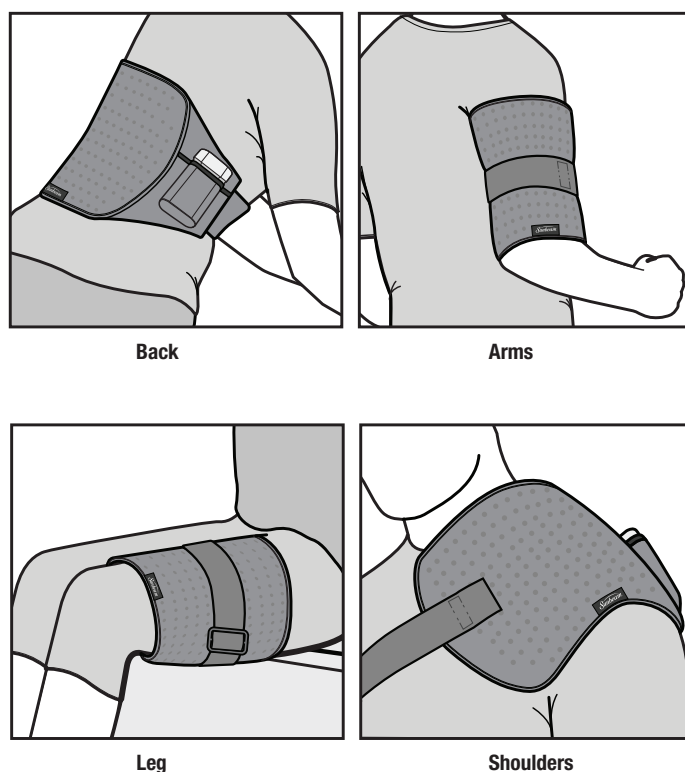HOW TO USE:

Charge It

- Before first use, fully charge the battery unit/heating pad controller by connecting charger cord into the charger port and plug into AC wall jack. Note: Charge in an ambient temperature between 41°F and 95°F (5°C and 35°C).
- Battery light will flash green during charging (See table below for indicators).
- Once battery light is solid green, the unit is fully charged. Charging may take up to 3 hours (using the standard charger) and up to 1 hour and 15 minutes (using the quick charger).

Select Heat Setting

- Use the power button to select heat settings by clicking consecutively to scroll through each setting. (See Figure 3 for button location)
- There are heat setting lights on battery unit to indicate if the unit is set to low, medium or high. (See Figure 4 for heat setting lights) Recommended treatment time is 20 minutes on high heat setting or 2 hours on low heat setting.

Apply Heating Pad

- Ensure the heating pad is connected to the battery unit by securely attaching the battery unit/connector to the pad connector before using.
- Put the battery unit into the pocket, controller side up.
- Place lighter colored fabric side (Figure 1, number 1) flat against the body.
- Take the edge of the adjustable strap and fasten directly to the darker, exterior side fabric, using the fastening tab (Figure 1, number 3) on the end of the strap. The fastening tab can attach to any part of the exterior side fabric. (See Figure 2 for the best way to wrap each body part) NOTE: Be careful not to overlap the heated area of the pad.
- Once in place, tighten around the body by adjusting the strap using the slider for a comfortable and snug fit.
- To remove pad, detach the fastening tab by pulling up on the end of the strap. Unwrap pad from body.

After Use

- Disconnect pad cord from the battery unit by unplugging cord at the bottom of the pocket.
- Place charger on the lighter fabric side, roll pad around charger, wrap strap around pad, and secure with the tab for neat and tidy storage.
- If needed, recharge battery unit prior to next use – only use the power supply that was provided with the unit.
- Store battery unit and product in cool, dry location between uses.

CAUTION:

- DO NOT SIT ON, LEAN AGAINST, CRUSH OR FOLD PAD DURING USE.
- DO NOT OVERLAP THE PAD.
- DO NOT USE WHILE SLEEPING.
- SEE ADDITIONAL WARNINGS IN THE IMPORTANT SAFETY INSTRUCTIONS SECTION.

BATTERY UNIT:

- The battery can last up to 4 hours on low heat setting and up to 2 hours on high heat setting.
- The battery light indicator will turn red or green when the battery is either low, charging, or fully charged. (See table 1 for the battery light colors per setting)
- Battery runtime will reduce over time. Battery life depends on usage and storage conditions.

Table 1:

Battery Light Indicator:	
Solid Red	Battery is low <20%
Flashing Red	Battery is very low <10%
Solid Green	Battery is fully charged
Flashing Green	Battery is plugged in and charging
Three Time Flash of Three White Dots	<ul style="list-style-type: none"> Battery has power, but pad is not connected to the battery unit Battery unit is attempted to be used while connected to the standard charger. Note: If using quick charger, pad can be used while charging.

HOW TO USE YOUR HEATING PAD WITH MOIST HEAT THERAPY

In some cases, you may find that moist heat therapy provides more comfort and relief. When using moist heat therapy, you should protect furniture and clothing from becoming damp.

- Dampen the lighter colored fabric side (Figure 1, number 1) of the heating pad with water by patting the surface with a wet sponge or cloth, or use a spray bottle to apply a fine mist. Do not thoroughly soak pad.
- Do not spray or dampen the pad connector or battery unit. (Figure 5)

DANGER: DO NOT IMMERSE THE PAD OR WIRING IN WATER OR PLACE UNDER RUNNING WATER. FAILURE TO COMPLY MAY RESULT IN PRODUCT FAILURE AND/OR PERSONAL INJURY. THE ELECTRICAL PAD CONNECTOR MUST BE THOROUGHLY DRY BEFORE ESTABLISHING ELECTRICAL CONNECTION.

Figure 3


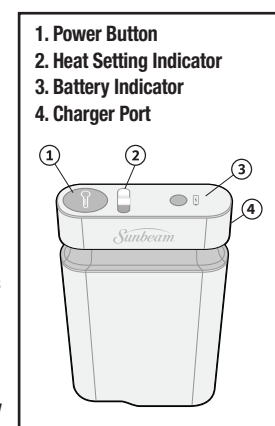
Figure 4

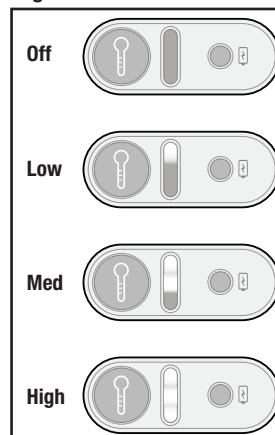WASHING INSTRUCTIONS

DO NOT DRY CLEAN! DO NOT USE BLEACH! DO NOT USE WRINGER! DO NOT IRON!

- Disconnect pad from charger if using quick charger while in use. DO NOT PULL BY THE CORD.
- Remove battery unit/heating pad controller from pocket and disconnect pad cable from unit.
- Set battery unit aside. DO NOT WASH.
- Tuck pad cable into pocket to protect it during the wash cycle.
- Machine Wash pad in cold water in gentle cycle or hand wash. Tumble dry on low heat and remove promptly from dryer or drape to air dry.
- DO NOT plug in pad connector to battery unit until pad and connector are completely dry.

For information on service or product questions, please visit www.sunbeam.com or contact customer service 1-800-435-1250

Scan for operating instructions, and warranty details.

How to scan QR Code:

- Open your camera app or QR Code scanning app.
 - Position your phone so the QR Code appears in the viewfinder.
 - Link will appear at the top.
- Or you can visit:
www.sunbeam.com/heatingpads/instructionmanuals

LIMITED WARRANTY

Sunbeam Products, Inc. ("Sunbeam") or if in Canada, Newell Brands Canada ULC ("NBC") warrants that for a period of one (1) year from the date of initial retail purchase, this product will be free from defects in material and workmanship, when used in accordance with the instructions provided.

Sunbeam (or NBC), at its option, will repair or replace this product or any component of the product found to be defective during the warranty period.

Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. Do NOT attempt to repair or adjust any electrical or mechanical functions on this product. Doing so will void this warranty.

This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. Sunbeam (or NBC) dealers, service centers, or retail stores selling Sunbeam (or NBC) products do not have the right to alter, modify or in any way change the terms and conditions of this warranty.

This warranty does not cover normal wear or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than Sunbeam (or NBC) or an authorized Sunbeam (or NBC) service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

What are the limits on Sunbeam's (or NBC's) Liability?

Sunbeam (or NBC) shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition. Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty. Sunbeam (or NBC) disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise. Sunbeam (or NBC) shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party. Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you. This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

How to Obtain Warranty Service

In the U.S.A.
If you have any question regarding this warranty or would like to obtain warranty service, please call 1 800 435-1250 and a convenient service center address will be provided to you.

In Canada
If you have any question regarding this warranty or would like to obtain warranty service, please call 1 800 667-8623 and a convenient service center address will be provided to you.

In the U.S.A. this warranty is offered by Sunbeam Products, Inc. located in Boca Raton, Florida 33431.

In Canada, this warranty is offered by Newell Brands Canada ULC, located at 208 Herford Street, Brampton, Ontario L6Y 0M1. If you have any other problem or claim in connection with this product, please write our Consumer Service Department.

PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE.

FCC CERTIFIED

The device complies with part 15 of the FCC rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This product has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This product generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures: Reorient or relocate the receiving antenna. Increase the separation between the product and receiver. Connect the product into an outlet on a circuit different from that to which the receiver is connected. Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

KEEP DRY
(Cleaning instructions: controller do not wet / must be kept dry)

RECYCLE
Controller Contains Lithium Battery

READ INSTRUCTIONS

DISPOSE OF PROPERLY
Not for general waste disposal

This product has earned the UL Listed Mark

Note: As of October 2019, this lithium battery-powered appliance meets FAA requirements for carry-on and checked baggage. Please check the FAA website for most up to date regulations.

©2020 Sunbeam Products, Inc. All rights reserved.

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